
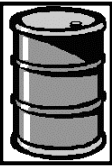




Product name:	MobileDiagnost wDR	
Identification code(s)	9890-010-89522	
Total weight (in Kg)	565 - 580 Kg	
Producer/ Manufacturer	Name company:	Philips Medical Systems DMC GmbH
	Address:	Roentgenstrasse 24
	Zip code:	22335 Hamburg
	Country:	Germany

Recycle Info	Items:	Location
Special attention 	Discharge Power Capacitors before handling the unit	Fig. 4
	Thallium surface on glass plate of X-Ray Detector	Fig. 3


	Items:	Location
 Fluids / Gases	Oil tank / Remove oil before dismantling	Fig. 1,4
	HV-transformer oil contains no PCBs	

	Type:	Location
Batteries  To be Removed	Lead Battery	Fig. 6
	Lithium Battery	Fig. 3 ; 5
	X-Ray Detector 3; Printed Circuit Board 5	

	Substances:	Location
Hazardous  To be Removed	Lead X-Ray tube 1; Collimator 2	Fig. 1 ; 2
	Thallium surface on glass plate of X-Ray Detector	Fig. 3

Show locations of materials mentioned on the previous sheet. (Picture information)

Mobile wDR with fix column m = 565 kg

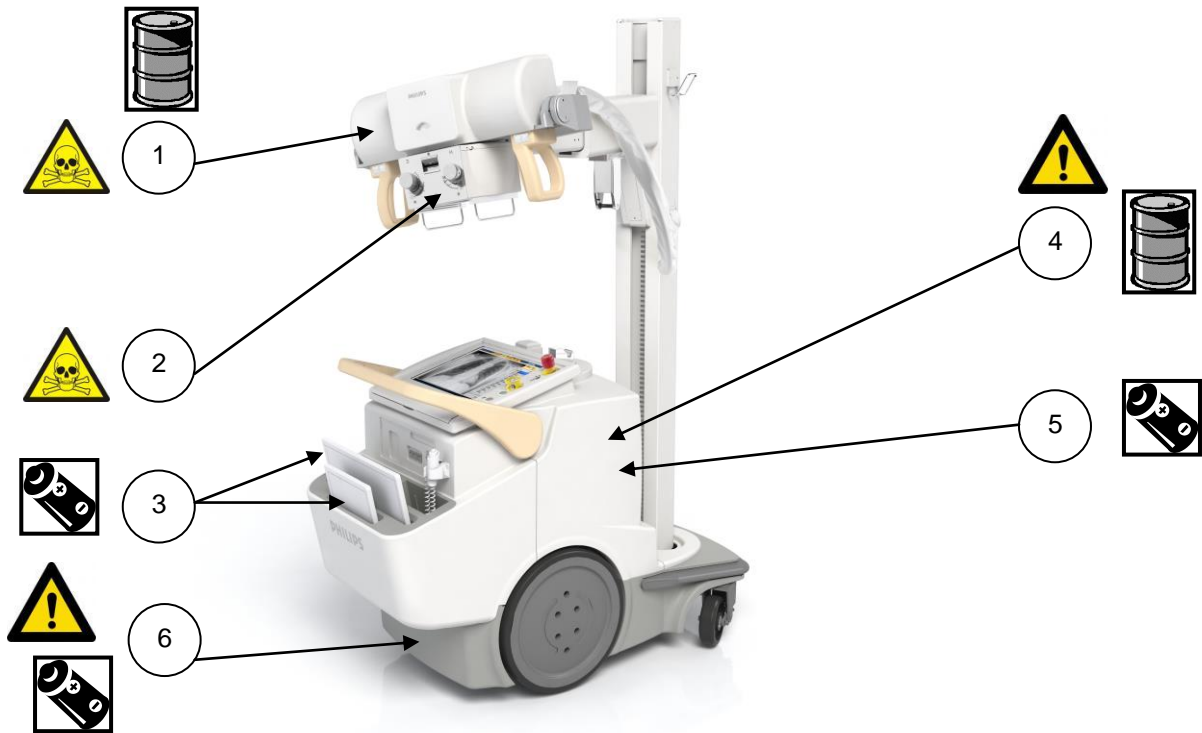

Mobile wDR with telescopic column m = 580kg



Recycling Passport Revision History

Rev	Description	Date
1.0	approved	23-Jun-2014
1.1	Add Mobile wDR 2.1 with telescopic column 580kg	24-Feb.-2016