
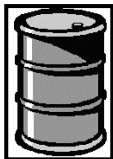




<b>PHILIPS HEALTHCARE</b>	<b>Product Recycling Passport DigitalDiagnost Rel. 4.x</b>	Page 1 of 4
-------------------------------	--	-------------

<b>Product name:</b>	<b>DigitalDiagnost Rel. 4.x</b>	
<b>Identification code(s)</b>	<b>9890-010-83904</b>	
<b>Total weight (in Kg)</b>	<b>1013 kg – 1515 kg</b>	
<b>Producer/ Manufacturer</b>	<b>Name company:</b>	<b>Philips Medical Systems DMC GmbH</b>
	<b>Address:</b>	<b>Röntgenstraße 24</b>
	<b>Zip code:</b>	<b>22335 Hamburg</b>
	<b>Country:</b>	<b>Germany</b>

<b>Recycle Info</b>	<b>Items:</b>	<b>Location</b>
<b>Special attention</b> 	<b>Lead in tube, Collimator, X-Ray Detector</b>	<b>Fig. 1 ; 2</b>
	<b>Thallium surface on glass plate of X-Ray Detector</b>	<b>Fig. 2 ; 5</b>
 <b>Fluids / Gases</b>	<b>Items:</b>	<b>Location</b>
	<b>X-ray tube oil tank / Remove oil before dismantling</b>	<b>Fig. 1</b>
	<b>X-ray generator oil tank not visible</b>	
<b>Batteries</b>  <b>To be Removed</b>	<b>Type:</b>	<b>Location</b>
	<b>Lithium Ion Battery in X-Ray Detector or Charger</b>	<b>Fig. 4</b>
	<b>Lithium Battery PC printed circuit board</b>	<b>Fig. 6</b>
	<b>Battery Remote Control</b>	<b>Fig. 7</b>
<b>Hazardous</b>  <b>To be Removed</b>	<b>Substances:</b>	<b>Location</b>
	<b>Lead in X-ray tube and Collimator</b>	<b>Fig. 1</b>
	<b>Lead in X-Ray Detector &amp; Grid</b>	<b>Fig. 5 ; 2</b>
	<b>Thallium surface on glass plate of X-Ray Detector</b>	<b>Fig. 5</b>

Ceiling Suspension CSM 9890-010-8870x



BuckyDiagnost VS1 9890-010-8364x



DigitalDiagnost TH2 9848-600-0222x

Vertical Stand 2 (VS2) 9897-010-0141x

Remote Control



Ceiling Suspension CSM 9890-010-8870x

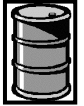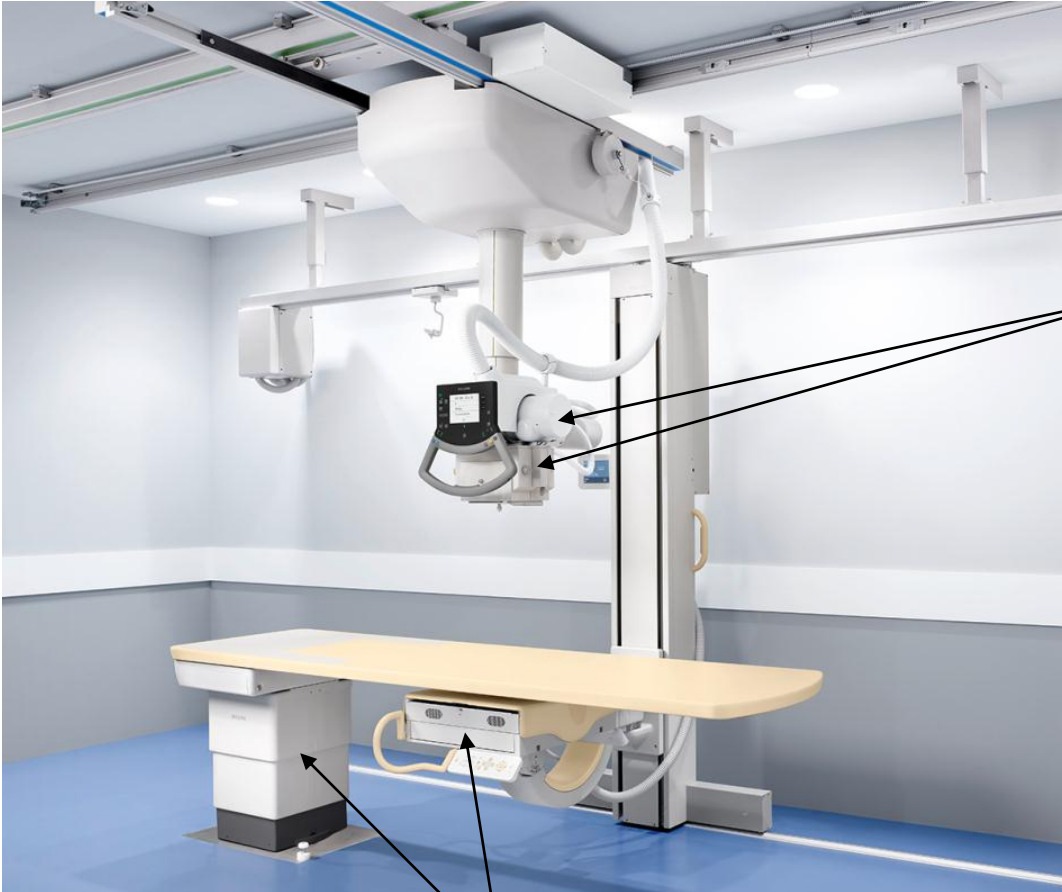


Single Sided Table TH-S



Option: SkyPlate 9897-010-00821

X-ray Detector SkyPlate


Charger for battery SkyPlate

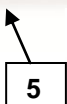

**Eleva Workspot Userstation for DigitalDiagnost  
9890-010-87551**



**Batteries Type lithium to be removed**

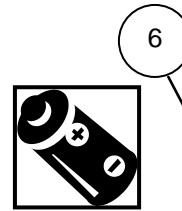

**Sub-System**

<b>X-Ray Tube assembly</b>	<b>9890-000-86481</b>	<b>see area X-Ray Tubes</b>
<b>X-Ray Generator M-Cabinet CXA</b>	<b>9890-010-88391</b>	<b>see area X- Generators</b>

**Recycling Passport Revision History**

<b>Ver.</b>	<b>Description</b>	<b>Date</b>
<b>00</b>	<b>DigitalDiagnost Rel. 4.0</b>	<b>18-Aug-2014</b>
<b>01</b>	<b>DigitalDiagnost Rel. 4.1</b>	<b>17-Nov-2014</b>