


PHILIPS HEALTHCARE	Product Recycling Passport Allura FD20/15	Page 1
-------------------------------	--	--------

Product name:	Allura Xper FD20/15	
Identification code(s)	722058	
Total weight (in Kg)	3000 Kg (approx)	
Producer /Manufacturer	Company:	Philips Healthcare
	Address:	Veenpluis 6
	Zip code:	5684 PC Best
	Country:	The Netherlands
	Electronic info:	DHF220780_Allura_FD20_15 recyclingPassport_V00_authorized.doc

Recycle info	Substances:	Location
Special attention 	Items:	Location
	Be aware of possible contaminated system parts and materials.	See Figure 1
	For dismantling activities Treatment Facilities must consider the national requirements.	
	For personnel that can come in contact with contaminated material, preventive measurements pursuant to national requirements must be taken into account.	
	Removal of units/weights can cause the system to tilt!	
	Removal of units/weights can cause unexpected movements.	
Hazardous  To be removed	Lead (Pb) for X-Ray Shielding	
	Beryllium copper (BeCu)	electronics
	Mercury (Hg) in specific LCD Screens	Monitors
	LEAD in X-Ray grid.	Grid
	Lead for soldering in Electronics, systems manufactured before July 2014.	All electronics.
Fluids / Gases 	Items: Tube cooling Oil. Shell Diala S3 ZX-IG Detector (when cooled) Glycoshell SF	Location X-Ray tube Detector

Note: To facilitate recycling, all plastic parts > 50 grams are marked according to ISO 11469 & ISO 1043

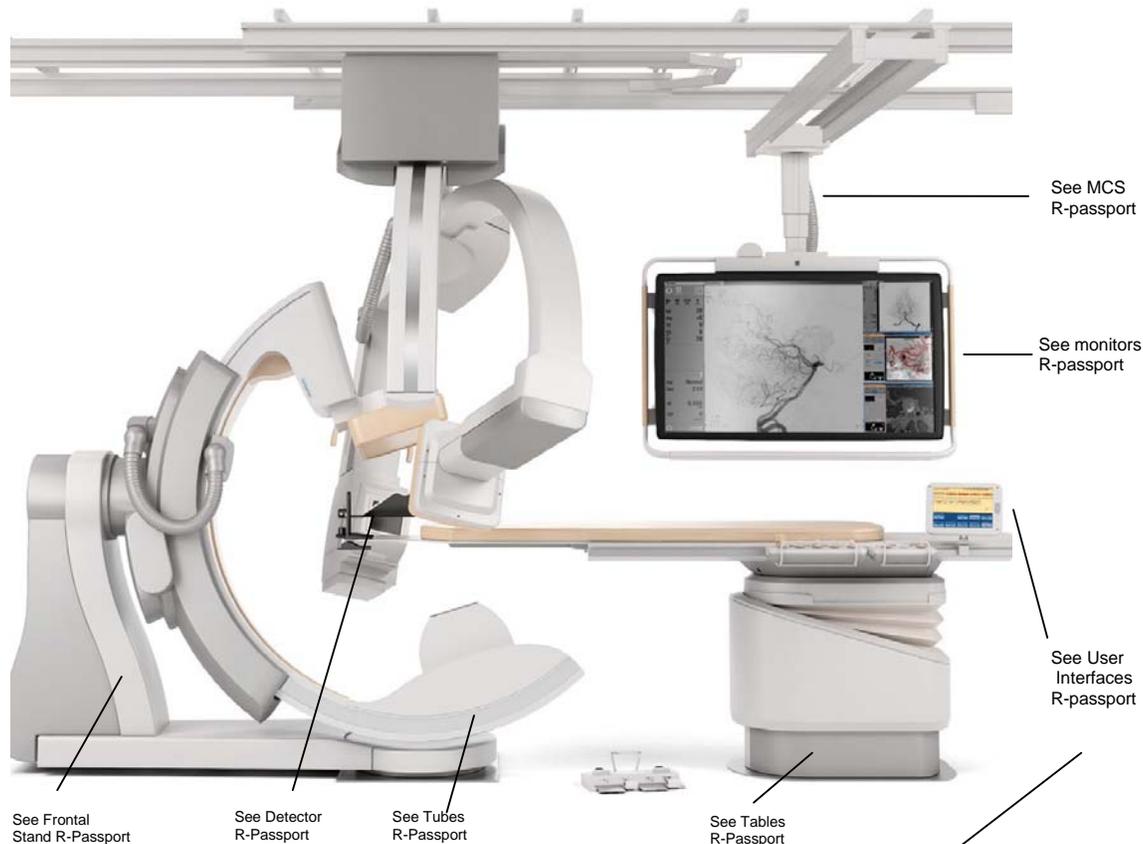

For Stand:

- See Frontal Stand recycling passport [includes FD20 configuration].
- See Lateral Stand Recycling passport. [Includes FD15 configuration]

Technical Room Cabinets:

- See M-Cabinet R-Passport
- See R-Cabinet R-Passport
- See X-Ray Generator Passport

Building blocks recycling passports are available on:
<http://www.philips.com/recycling>