

<b>Product name</b>	Recycling passport Mammotrak trolley 1.5T and 3.0T
Identification code(s)	98960301987x (1.5T version) and 98960302084x (3.0T version) (where x can be any digit from 1....9)
Total weight [kg]	142 kg
<b>Manufacturer</b>	<b>Philips Healthcare Nederland B.V.</b>
Address	Veenpluis 4-6 Best
Zip code	5684PC
Country	The Netherlands
Electronic info	APT591-08-4701 Mammotrak excl RF coils WEEE passport

## Recycling information

Special Attention	Items	Location
	Under the carrier in the spacer there is a cable slab. An RF cable, a nurse call cable, and the headset tube run through this cable slab.	Nr 1: Figure 2
	An ODU connector (for the RF cable) and connectors for the nurse call and headset are positioned under the lid at the head-end of the carrier	Nr 2: Fig 1 and fig 4
	A cable plus ODU connector and connection to the headset and nurse call are at the side of the spacer.	Nr 3: Fig 3 and 4 ODU Connector: fig 5.

Batteries	Items, type	Location
	Not applicable	

Hazardous substances	Items	Location
	Beryllium copper contact pins on ODU connectors	Fig 5

Name: Recycling Passport Mammotrak trolley 1.5T and 3T		12NC: 98960301987x (1.5T) and 98960302084x (3.0T)		Version: <b>00</b> Doc.Date: <b>2008-10-15</b>	
Author:		Pag.: <b>1 of 4</b>	Group No: 5790	Section: -	Status: <b>Authorized</b>

Fluids and gases	Items	Location
	Not applicable	

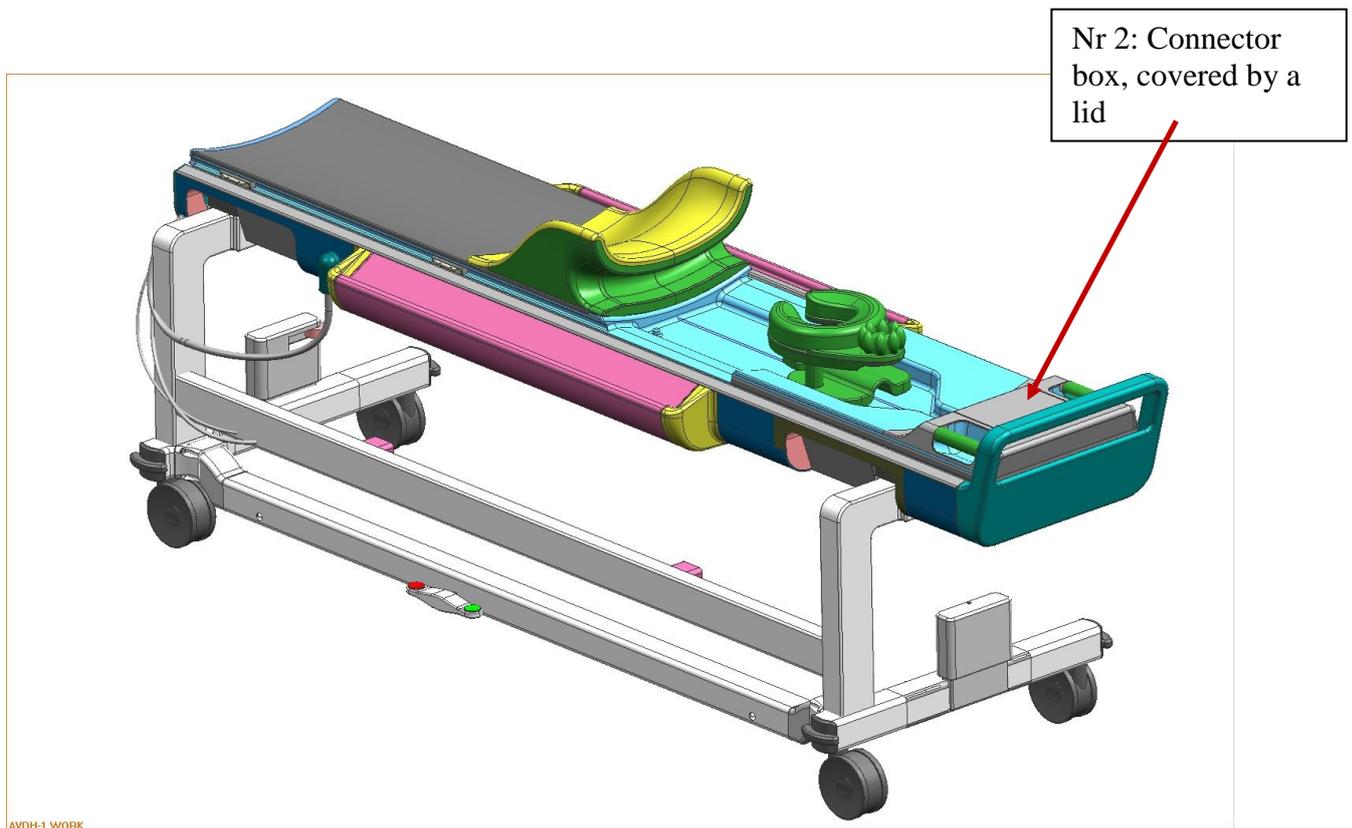


Fig 1: Drawing of Mammotrak with carrier, Ramp and Headrest

Name: Recycling Passport Mammotrak trolley 1.5T and 3T		12NC: 98960301987x (1.5T) and 98960302084x (3.0T)		Version: <b>00</b> Doc.Date: <b>2008-10-15</b>	
Author:		Pag.: <b>2 of 4</b>	Group No: 5790	Section: -	Status: <b>Authorized</b>
© PHILIPS HEALTHCARE NEDERLAND BV, BEST 2008			Ref. XJV032011.07 Rev.03		

Nr 1: Cable slab, containing RF cable, nurse call cable and tube for headset


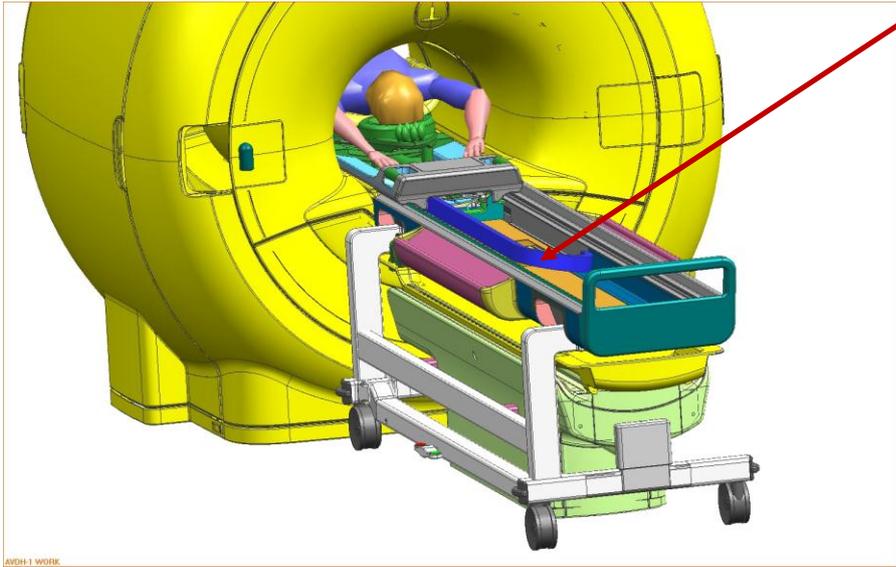


Fig 2: Artist Impression of Mammotrak docked to the MR scanner

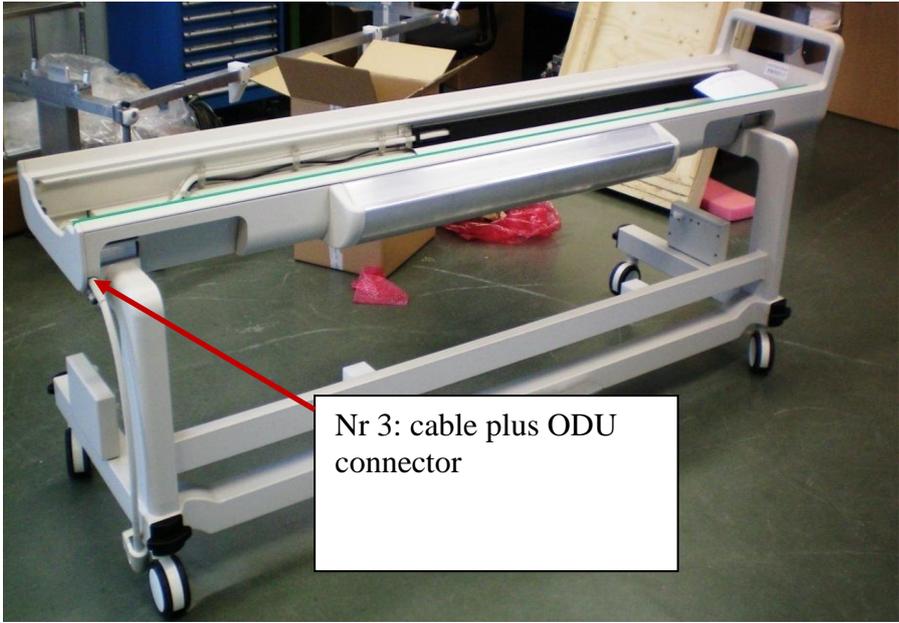

Nr 3: cable plus ODU connector

Fig 3: Mammotrak without carrier

Name: Recycling Passport Mammotrak trolley 1.5T and 3T		12NC: 98960301987x (1.5T) and 98960302084x (3.0T)		Version: <b>00</b> Doc.Date: <b>2008-10-15</b>	
Author:		Pag.: <b>3 of 4</b>	Group No: 5790	Section: -	Status: <b>Authorized</b>
© PHILIPS HEALTHCARE NEDERLAND BV, BEST 2008			Ref. XJV032011.07 Rev.03		

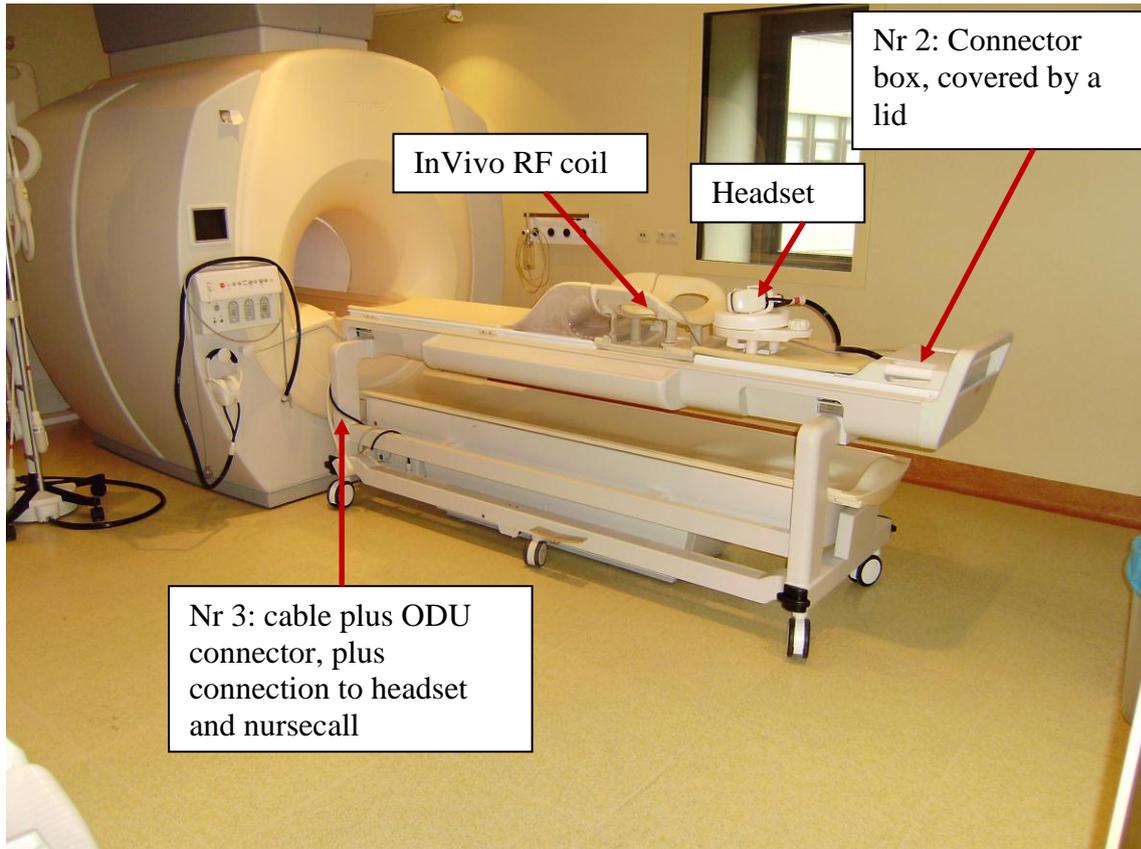


Fig 4: Mammotrak with headset and Invivo RF coil

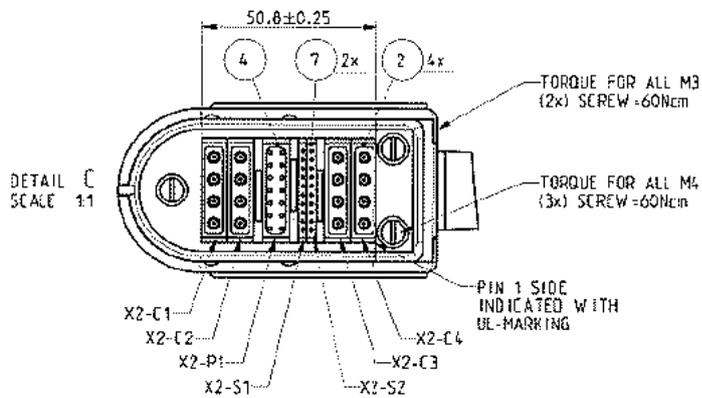


Fig 5: Layout of ODU connector

Name: Recycling Passport Mammotrak trolley 1.5T and 3T		12NC: 98960301987x (1.5T) and 98960302084x (3.0T)		Version: <b>00</b> Doc.Date: <b>2008-10-15</b>	
Author:		Pag.: <b>4 of 4</b>	Group No: 5790	Section: -	Status: <b>Authorized</b>
© PHILIPS HEALTHCARE NEDERLAND BV, BEST 2008			Ref. XJV032011.07 Rev.03		