

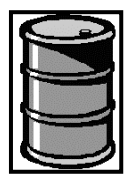



Product name:	Frontal Stand FD20 Ceiling / Floor version		
Identification code(s)	4522 901 00470		
Total weight (in Kg)	1000 Kg (approx)		
Producer /Manufacturer	Name company:	Philips Medical Systems	
	Address:	Veenpluis 6	
	Zip code:	5684 PC Best	
	Country:	The Netherlands	
	Electronic info:	Frontal Stand_FD20_579.doc	


Recycle info	Substances:	Location
Special attention 	Items:	Location
	Be aware of possible contaminated system parts and materials.	See Figure 1
	For dismantling activities Treatment Facilities must consider the national requirements.	
	For personnel that can come in contact with contaminated material, preventive measurements pursuant to national requirements must taken into account.	
	Removal of units/weights can cause the system to tilt!	
	Removal of units/weights can cause unexpected movements.	
	High voltage parts are marked with 	
 Fluids / Gases	Items:	Location
	Flat Detector; cooling fluid (glycol) may leak from tube and cooling hoses during and after dismantling.	Figure 1, item C
	X-ray tube; cooling fluid (oil) may leak from tube and cooling hoses during and after dismantling.	Figure 1, item D
Hazardous  To be removed	Lead (Pb) for soldering	All electronics
	Beryllium copper (BeCu)	Some electronics

Name:	WEEE Frontal Stand Allura XPER	12NC:	4522 901 00470	Version: 00
Author: W.Bleyen		Pag.: 1 of 4	Group No: 579	Doc.Date: 2006-09-01
			Section: -	Status: Approved


Note: To facilitate recycling, all plastic parts > 50 grams are marked according to ISO 11469 & ISO 1043

Frontal stand
Ceiling version

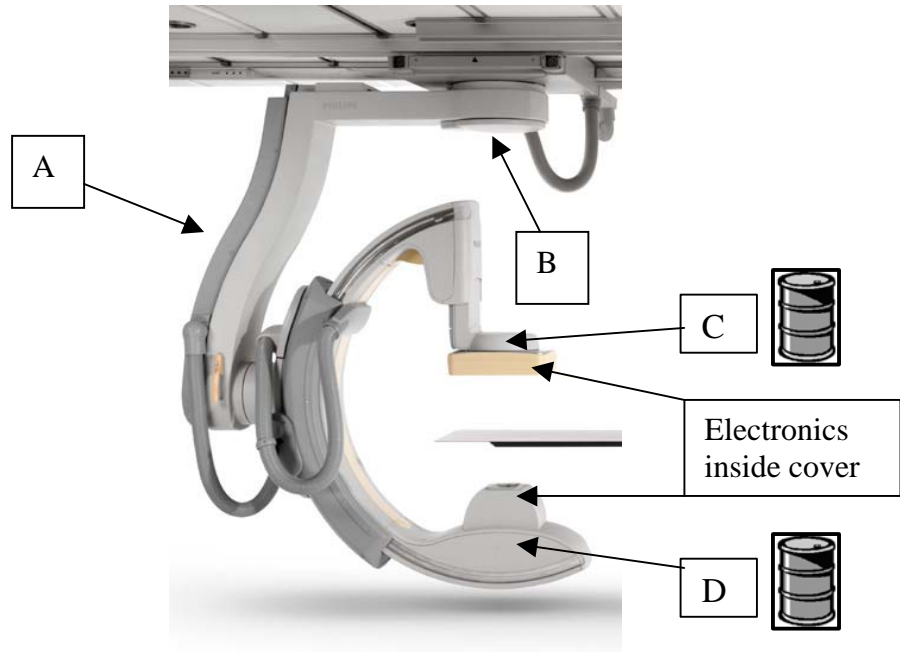Figure 1



View A

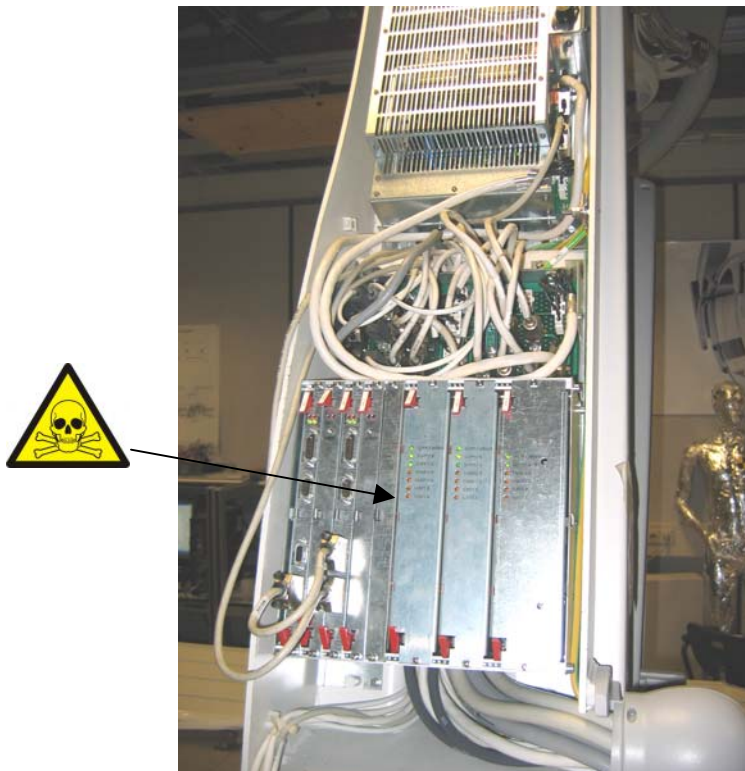


Name: WEEE Frontal Stand Allura XPER		12NC: 4522 901 00470		Version: 00 Doc.Date: 2006-09-01	
Author: W.Bleyen	Pag.: 2 of 4	Group No: 579	Section: -	Status: Approved	

View B




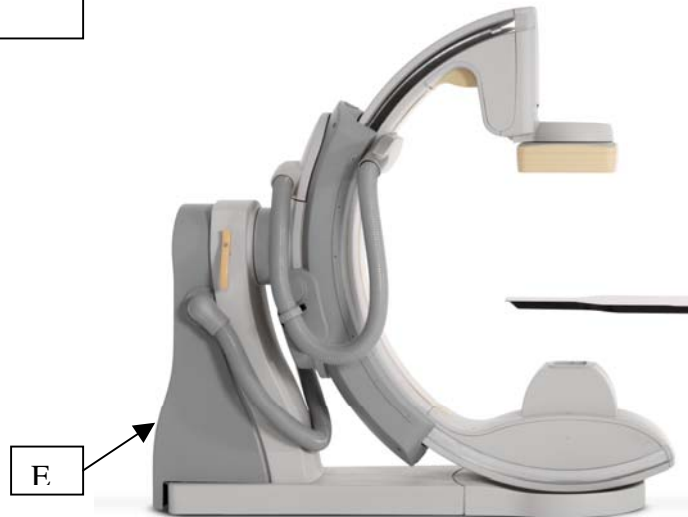Frontal stand
Floor version



Name:	WEEE Frontal Stand Allura XPER	12NC:	4522 901 00470	Version:	00
Author:	W.Bleyen	Pag.:	3 of 4	Doc.Date:	2006-09-01
		Group No:	579	Section:	-
				Status:	Approved

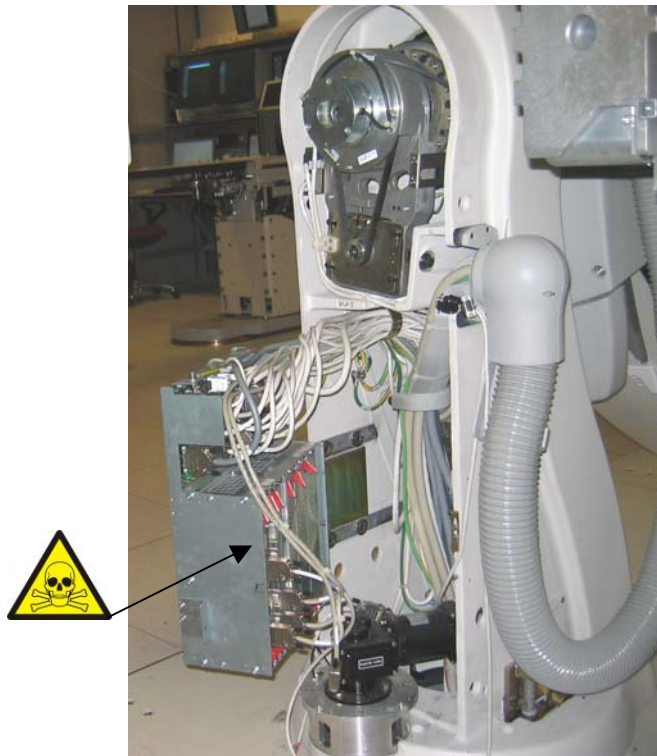View C



Recycling Passport Revision History

Rev	Description	Date
1	First Passport	Sept. 1, 2006

Recycling Passport Approval Signatures

Written By PMS Tom Verzaal
--

Reviewed By (optional)

Approved By PMS W.Bleyen
--

Approved By

Name: WEEE Frontal Stand Allura XPER		12NC: 4522 901 00470		Version: 00 Doc.Date: 2006-09-01	
Author: W.Bleyen	Pag.: 4 of 4	Group No: 579	Section: -	Status: Approved	