




<b>Product name:</b>	<b>PATIENT TABLE NL</b> <b>PATIENT TABLE NL/3T</b> <b>PATIENT TABLE I.V.</b> <b>PATIENT TABLE 1.0T/1.5 T WE</b> <b>PATIENT TABLE 3.0T WE</b> <b>PATIENT TABLE 1.5T I/T WE</b>		
<b>Identification code(s)</b>	<b>4522 131 6150 x</b> <b>4522 131 6160 x</b> <b>4522 131 9756 x</b> <b>4522 132 2100 x</b> <b>4522 132 2145 x</b> <b>4522 132 2146 x</b>		
<b>Total weight (in Kg)</b>	<b>185</b>		
<b>Producer /Manufacturer</b>	Name company:	Philips Medical Systems Nederland B.V.	
	Address:	Veenpluis 4-6	
	Zip code:	5680 DA BEST	
	Country:	The Netherlands	
	Electronic info:	PatientTable_452213161500_579_v01.doc	

<b>Recycle info</b>	<b>Substances:</b>	<b>Location</b>
<b>Special attention</b> 	<b>Items:</b>	<b>Location</b>
	These products contain heavy construction elements. Appropriate care shall be taken during disassembly.	1
	These products contain a scissor mechanism. Care shall be taken during disassembly to prevent injury due to collaps of the mechanism.	2
	The power supplies contain electrolytic capacitors. The power supplies contain Printed Circuit Boards.	3
	These products contain Printed Circuit Boards inside the indicated housings.	4
	These products have external and internal electric cabling.	5
	The actuator contains an electrolytic capacitor.	2
<b>Hazardous</b> 	The EMC housing inside these products, contains sticky fingers. The sticky fingers consist of beryllium copper, a hazardous material.	6
	The Printed Circuit Boards contain lead based solder.	3, 4
<b>To be removed</b>		

Name: <b>Recycling Passport Patient Tables</b>		12NC: <b>4522 131 6150 0</b>	Version: <b>00</b> Doc.Date: <b>2006-09-14</b>
Author: <b>G.J. Lenting</b>	Pag.: <b>1 of 3</b>	Group No: <b>5790</b>	Section: <b>-</b> Status: <b>Authorized</b>

Locations as mentioned in the Passport ( Pictures information)

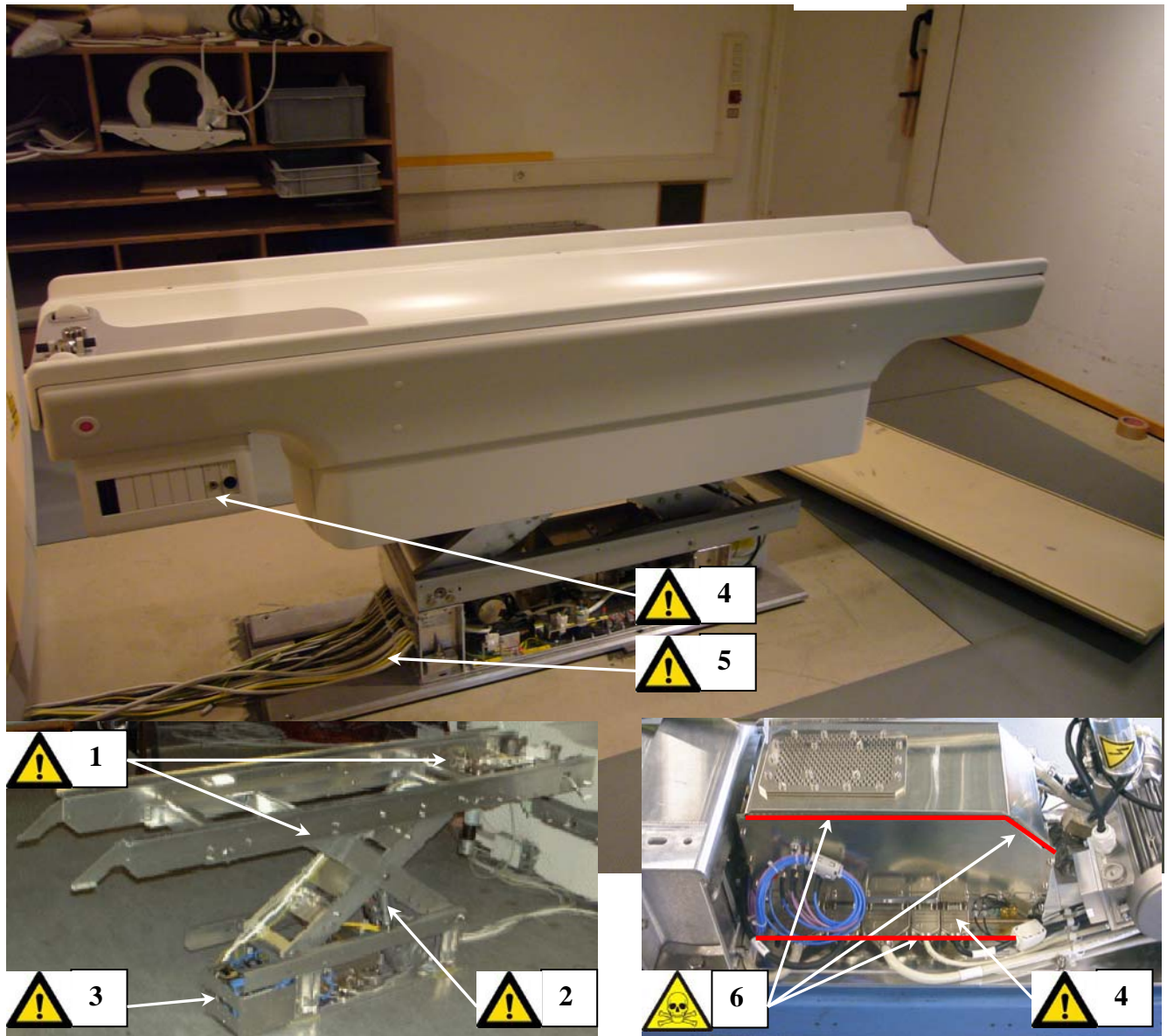

Name: <b>Recycling Passport Patient Tables</b>		12NC: <b>4522 131 6150 0</b>	Version: <b>00</b> Doc.Date: <b>2006-09-14</b>
Author: <b>G.J. Lenting</b>	Pag.: <b>2 of 3</b>	Group No: <b>5790</b>	Section: <b>-</b> Status: <b>Authorized</b>

### Recycling Passport Revision History

Rev	Description	Date
00	First Passport	September 14 <sup>th</sup> , 2006

### Recycling Passport Approval Signatures

Written By GJ Lenting  
**GJ Lenting**

Reviewed By T Pessers  
**T Pessers**

Approved By

Approved By W Evers  
**W Evers**

Name: <b>Recycling Passport Patient Tables</b>			12NC: <b>4522 131 6150 0</b>			Version: <b>00</b> Doc.Date: <b>2006-09-14</b>			
Author: <b>G.J. Lenting</b>		Pag.: <b>3 of 3</b>		Group No: <b>5790</b>		Section: <b>-</b>		Status: <b>Authorized</b>	