






**PHILIPS MEDICAL SYSTEMS / PRODUCT RECYCLING PASSPORT**

Printed copies are uncontrolled

<b>Product name:</b>	<b>Precedence (SPECT)</b>	
<b>Identification code(s)</b>	2169-3000A, 2169-3001A, 2169-3002A, 2169-3003A	
<b>Total weight</b>	<b>9000 lbs (approximately)</b>	
<b>Producer/ Manufacturer</b>	Name company:	ADAC Laboratories, A Philips Medical Systems Company
	Address:	540 Alder Dr., Milpitas, CA
	Zip code:	95035
	Country:	United States

<b>Recycle Info</b>	<b>Items:</b>	<b>Location/(Qty)</b>
<b>Special attention</b> 	Be aware of possible contaminated system parts and materials.	
	For dismantling activities Treatment Facilities must consider the national requirements.	
	For personnel that can come in contact with contaminated material, preventive measurements pursuant to national requirements must be taken into account.	
	Removal of units/weights may cause the system to tilt.	
	Removal of units/weights can cause unexpected motion.	
	High voltage parts are marked with 	
<b>Fluids / Gases</b>	<b>Items:</b>	<b>Location</b>
	None	
<b>Batteries</b>	<b>Type:</b>	<b>Location</b>
 <b>To be Removed</b>	Lithium Batteries on PC Boards and LCD Display	Figures 1 & 3
	Batteries for UPS	Figure 1
	Batteries for Hand controller	Figure 1


**PHILIPS MEDICAL SYSTEMS / PRODUCT RECYCLING PASSPORT**

<b>Hazardous</b>	<b>Substances:</b>	<b>Location/(Qty)</b>
	Lead (Pb) – Detector & Collimator	Figure 2
	Lead (Pb) for soldering	Figures 1& 3 All electronics (including PCBs, power supplies)
	Sodium Iodide & Thallium	Figure 2
	Lead (Pb) – Table PCBs	Figure 2

**NOTE: Recycle information for CT Gantry is covered under CT Brilliance Recycle Passport, No. 4535 671 49361**

Printed copies are uncontrolled

Printed copies are uncontrolled

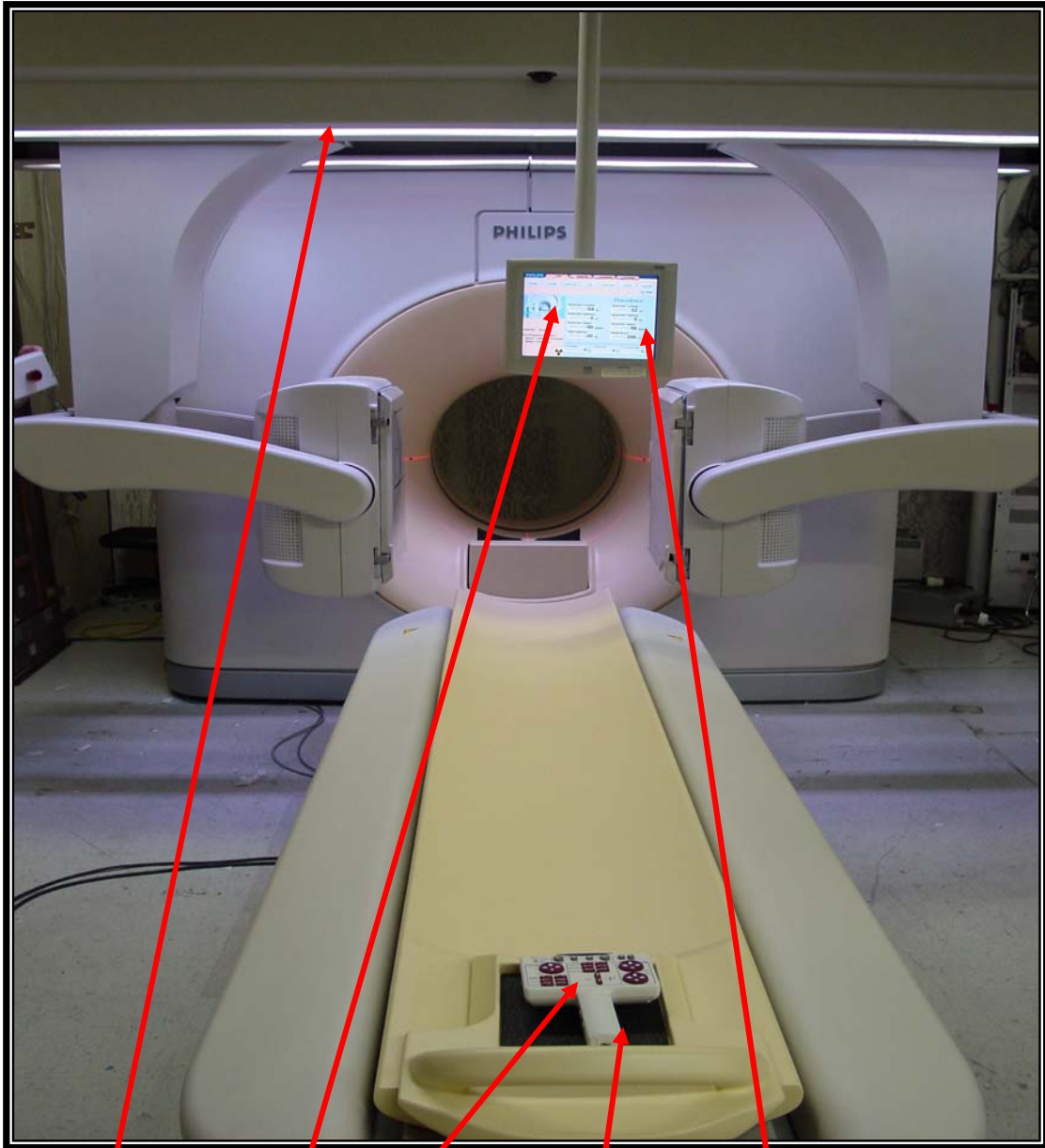

**Figure 1**  
Overall System



Batteries for UPS



Pb PCBs  
(inside LCD Display & Hand  
Controller)




Batteries for Hand Controller



Lithium Battery inside LCD Display

**Title: Precedence (SPECT) Recycling Passport**  
**Recycling Passport Number: 9756-0070 Rev: #1**

Printed copies are uncontrolled

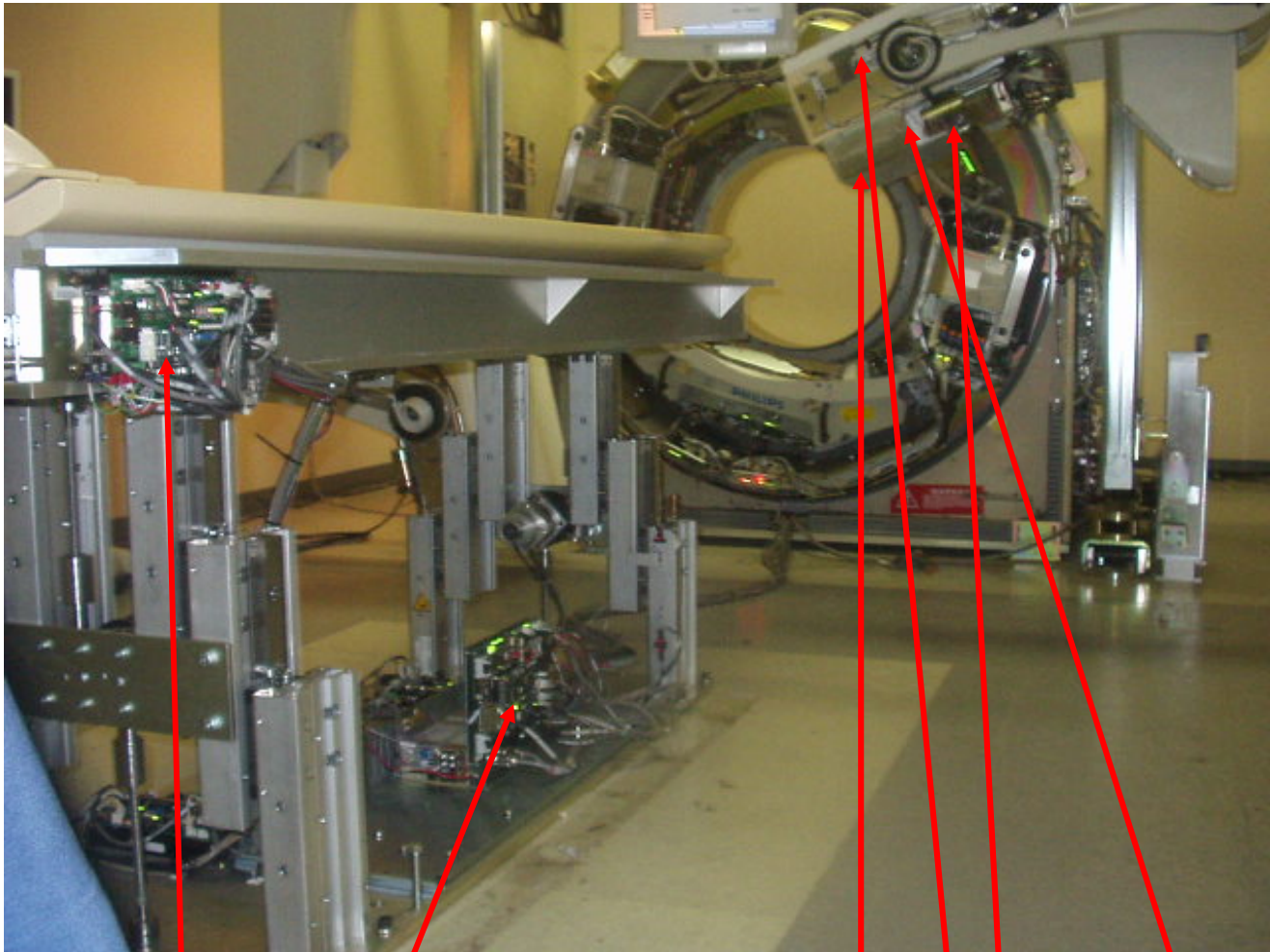

**Figure 2**  
**Internal View**



Pb  
PCBs



Sodium Iodide & Thallium  
(inside Detector Crystals)




Pb PCBs



Pb  
Detector &  
Collimator

Printed copies are uncontrolled

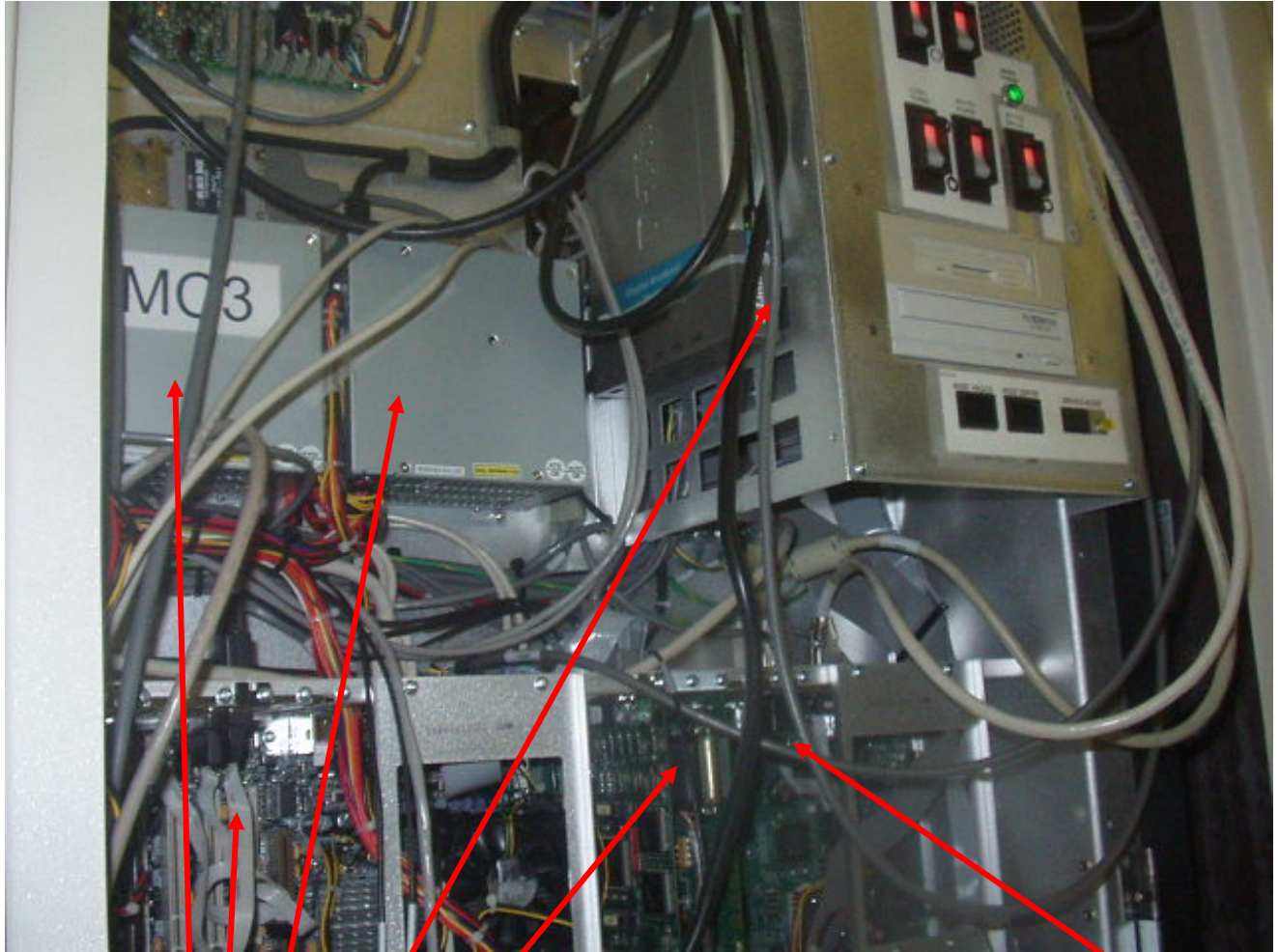

**Figure 3**  
**Internal View of SPECT Electronics Cabinet**



Pb  
PCBs, power  
supplies



Lithium Battery on  
PC Board