





Product name:	7T Achieva MR Scanner FCC Cabinet	
Identification code	9896 030 5052*	
Total weight (in Kg)	Approximately 256 kg	
Producer /Manufacturer	Name company:	Philips Medical Systems Cleveland
	Address:	595 Miner Rd., Highland Hts, Ohio
	Zip code:	44143
	Country:	USA
	Electronic info:	http://www.medical.philips.com/main/company/sustainability/recycling/mri/

Recycle info	Substances:	Location
 Hazardous To be removed	Lead (Pb) (solder – PC boards)	Figures 1, 2
	Beryllium-Copper EMI shielding gasket.	Figure 2
 Batteries To be removed	Type:	Location
	Li 3V button cell battery	Figure 1
	Magnet EDU: EDU – E7002 lead acid; EDU monitor 7107 NiCad	Figure 1
 Special attention	Item	Location
	Compressed air	Figure 2
	Storage capacitors.	Figure 2
	Magnet Monitoring Unit is not recyclable. Return to: Magnex Scientific Limited 6 Mead Rd, Oxford Ind. Park Yarnton, Oxford, UK OX5 1QU	Figure 1
 Fluids / Gases	Item	Location
	Compressed air line	Figure 2

* Denotes revision level wildcard designation. (Numbers range from 1-9).

Material Content (approximate)	Weight (Kg)	Material content (approximate)	Weight (Kg)
Iron (Fe)	154.1	Gold (Au)	-
Aluminum (Al)	42.4	Tungsten (W)	-
Copper (Cu)	23.3	Other	4.7
Plastics Epoxy etc	4.7		
Printed Circuit Boards	27.1		

Locations as mentioned in the Passport (Pictures information)

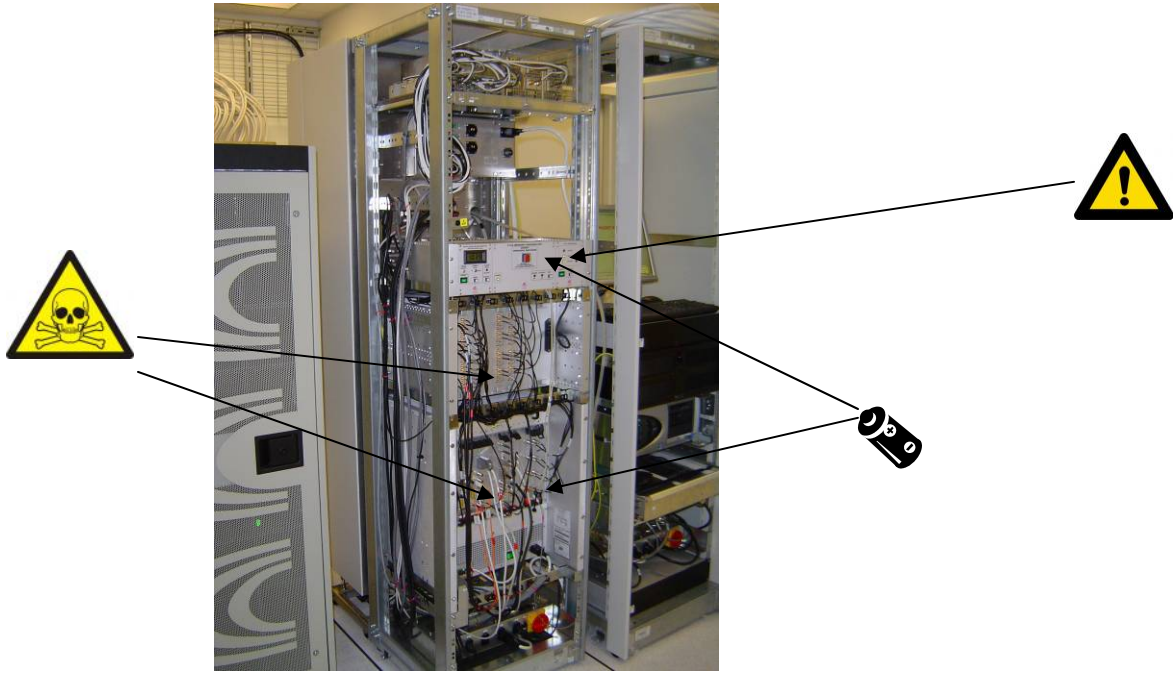


Figure 1. 7T FCC cabinet, front view.

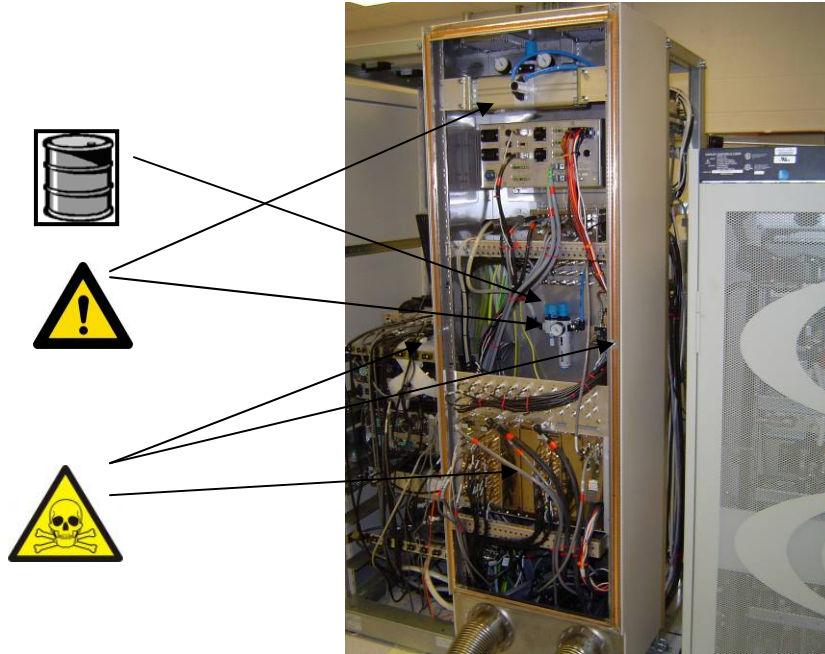


Figure 2. 7T FCC Cabinet, rear view.