



Log on to a better way to work

EncoreAnywhere patient management system

PHILIPS

RESPIRONICS

Discover how EncoreAnywhere can help you:

Save time

- Patient information is automatically downloaded
- No more faxing reports
- No more chasing down Encore SmartCards

Coordinate patient care more easily and more timely

- Share information with clinicians, physicians and sleep testing personnel via the secure Web site
- Implement prescription changes quickly and easily over the Internet

Offer your customers the latest advancement in sleep therapy

- Patients appreciate automatic downloads and quicker prescription changes
- Referrals appreciate ease-of-use, office efficiencies, and improved patient care

Demonstrate the unique benefits to your referrals

- Get sales support with our new Take Ten marketing program. In just ten minutes, your sales team can demonstrate the benefits of EncoreAnywhere to your referral sources.
- Get training support with our Take Ten setup packet that provides all of the tools needed to help your referrals learn to use the EncoreAnywhere system

See at a glance how EncoreAnywhere works:



User profiles determine icon notification availability on the My Day page.

The **My Day** page is a user's first stop to review patient treatment status.

The My Day page incorporates icons that automatically identify problem patients and their specific issues. With this information, users can quickly and easily prioritize patient intervention needs.



This icon indicates efficacy issues.



This icon indicates patient use issues.



This icon indicates user communication to care team members.



This icon indicates a prescription change has been initiated.



This icon indicates device communication issues.

Through EncoreAnywhere, the entire care team – including physicians and clinicians – has access to the same information and can communicate via notification messages. A comprehensive, collaborative approach to patient care helps ensure patient compliance and positive treatment outcomes.





Click through to a patient record to evaluate and respond to a patient's treatment needs in less time with more information.

The **Patient Summary** page provides patient information and a quick view of therapy status. Patient data is downloaded automatically through a wired modem or it is uploaded from an Encore SmartCard to the EncoreAnywhere system.

On the **Prescription** page the prescription can be viewed and changed. The prescription change will be uploaded automatically to the device through the wired modem within a preset time frame or can be sent via an Encore SmartCard.

Therapy data is easily viewed on the **Therapy Data** page by a specific time frame. Summary and daily detail reports can be easily generated.

The **Notes** page and alerts enable improved communication and collaborative patient care among providers, physicians and sleep testing facilities.

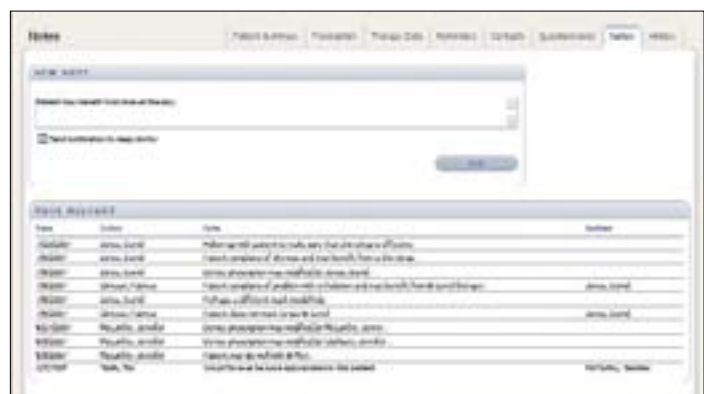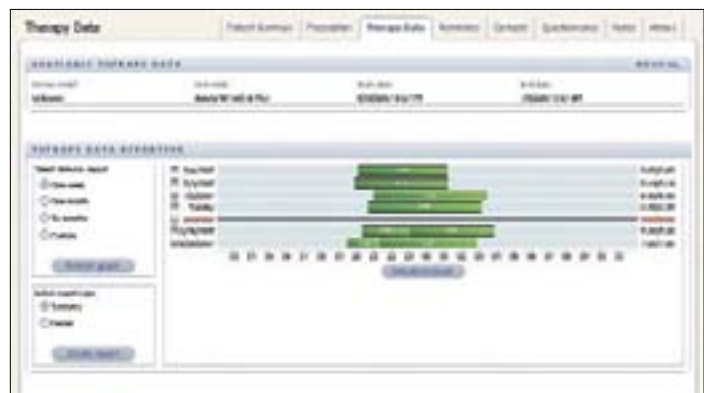

Real-time data access

The new M Series wired modem module enables real time access to patient data. The module fits into the standard M Series accessory slot in place of the SmartCard module. It will store patient data and download automatically to EncoreAnywhere based on a preset call program. The toll-free communication connection will not interfere with a patient phone line and will automatically upload pressure prescription changes.

The M Series wired modem module includes the following:

- M Series wired modem module
- Phone line/cord (14 feet)
- Phone line splitter
- One year of data downloads on a pre-programmed schedule



Extreme security

EncoreAnywhere is rigorously designed using the industry's best practices for data management. Philips Respironics has obtained an independent review and validation of the security measures protecting your data. Our customer service or sales representatives can provide you with written confirmation of the independent review and validation at your request.

Ordering information	
EncoreAnywhere	Part #
EncoreAnywhere subscription	1999999
EncoreAnywhere wired modem module	1037688
M Series wired modem call plan renewal	1999998
Take Ten referral marketing package	1056397
Take Ten referral setup and training package	1056398
Accessories	Part #
Encore USB SmartCard reader/writer	1014352
Encore SmartCard PCMCIA reader	1014353
Pack of Encore SmartCards (10 cards)	1003304
FOSQ clinical reprint	1000141
FOSQ questionnaire, English	1004550
Sleep therapy wired modem systems	Part #
REMstar Pro w/C-Flex and wired modem (includes humidifier)	DS400HW
REMstar Auto w/A-Flex and wired modem (includes humidifier)	DS510HW
BIPAP Auto w/Bi-Flex and wired modem (includes humidifier)	DS700HW

System requirements

Modem:

- Compatible with all M Series sleep therapy devices capable of using a SmartCard accessory module
- Land-based phone line (not wireless) with working phone service
- Phone jack located within 14 feet of the therapy device

Computer:

- Minimum 4 CPU and 512MB RAM
- Internet connectivity (not dial-up)
- Internet browser (IE 5.0 or higher)
- Adobe Reader (7.0 or higher)

Good for you and your customers

EncoreAnywhere is a new Web-based patient management system from Philips Respironics that's packed with benefits for you, your patients and your referrals.

Successful outcomes for sleep therapy patients can best be achieved through the collaborative effort of a team of healthcare providers. And, for the team to quickly respond to patient needs, it must have easy access to timely and accurate compliance information. With EncoreAnywhere, you and the entire healthcare team can view timely, accurate patient information, change prescriptions, and communicate easily – all through a secure Web site that's accessible from any computer, anywhere.

And the good news is it's easy to use.

With our new, wired modem modules, designed to work with Philips Respironics M Series sleep therapy devices, EncoreAnywhere automatically downloads data for instant access. What could be simpler?

Even better, EncoreAnywhere takes just minutes to learn.

In fact, in as little as ten minutes, your referrals can be up and running – saving time, saving money, and improving patient care.



**Philips Healthcare is part of
Royal Philips Electronics**

How to reach us

www.philips.com/healthcare

Fax: +31 40 27 64 887

Asia

+852 2821 5888

Europe, Middle East, Africa

+31 40 27 87246

Latin America

+55 11 2125 0764

North America

+1 425 487 7000

+1 800 285 5585 (toll free, US only)

Philips Respironics
1010 Murry Ridge Lane
Murrysville, PA 15668

Customer Service

+1 724 387 4000

+1 800 345 6443 (toll free)

Respironics Europe, Middle East, Africa

+33 1 47 52 30 00

Respironics Asia Pacific

+852 3194 2280

www.philips.com/respironics

Respironics, A-Flex, Bi-Flex, BiPAP, Encore, EncoreAnywhere, REMstar and SmartCard are trademarks of Respironics, Inc. and its affiliates. The C-Flex mark is used under license.



©2012 Koninklijke Philips Electronics N.V. All rights are reserved.

Philips Healthcare reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

CAUTION: US federal law restricts this device to sale by or on the order of a physician.

edoc SB 07/03/12 MCI 4104995