





PHILIPS HEALTHCARE	Product Recycling Passport	Page 1 of 7
-------------------------------	-----------------------------------	-------------

Product name:	Brilliance iCT Gantry	
Identification code(s)	4535 674 49901, 4535 675 31111	
Total weight (in Kg)	2520	
Producer/ Manufacturer	Name company:	Philips Medical Systems
	Address:	595 Miner Road Highland Heights, Ohio
	Zip code:	44143
	Country:	United States of America

Recycle Info	Items:	Location
Special attention 	Spring loaded cover	Figure 1 (1)
	Rotor unbalanced when a component is removed.	Figure 2 (2)
	Hazardous Voltage Warning	Figure 3 (9)
Fluids / Gases	Items:	Location
	Hydrocarbon Oil	Figure 2 (3)
	Pressure Springs: Must drill 2mm hole to relieve pressure before disposing of springs.	Figure 2 (4)
Batteries	Type:	Location
		
To be Removed		
Hazardous	Substances:	Location
	Printed Circuit Boards	Figure 2 (5) (11), Figure 2 (6), Figure 4 (10), Figure 5 (11), Figure 2 (8), Figure 2 (7)
	Collimator-Lead	Figure 2 (7)
To be Removed		

Title: Recycle Passport Brilliance iCT Gantry
 Recycling Passport Number: 453567534481 Rev: B

Printed copies are uncontrolled

Show locations of materials mentioned on the previous sheet. (Picture information)

Printed copies are uncontrolled



(1)



Figure 1: Front Cover of Gantry

Printed copies are uncontrolled


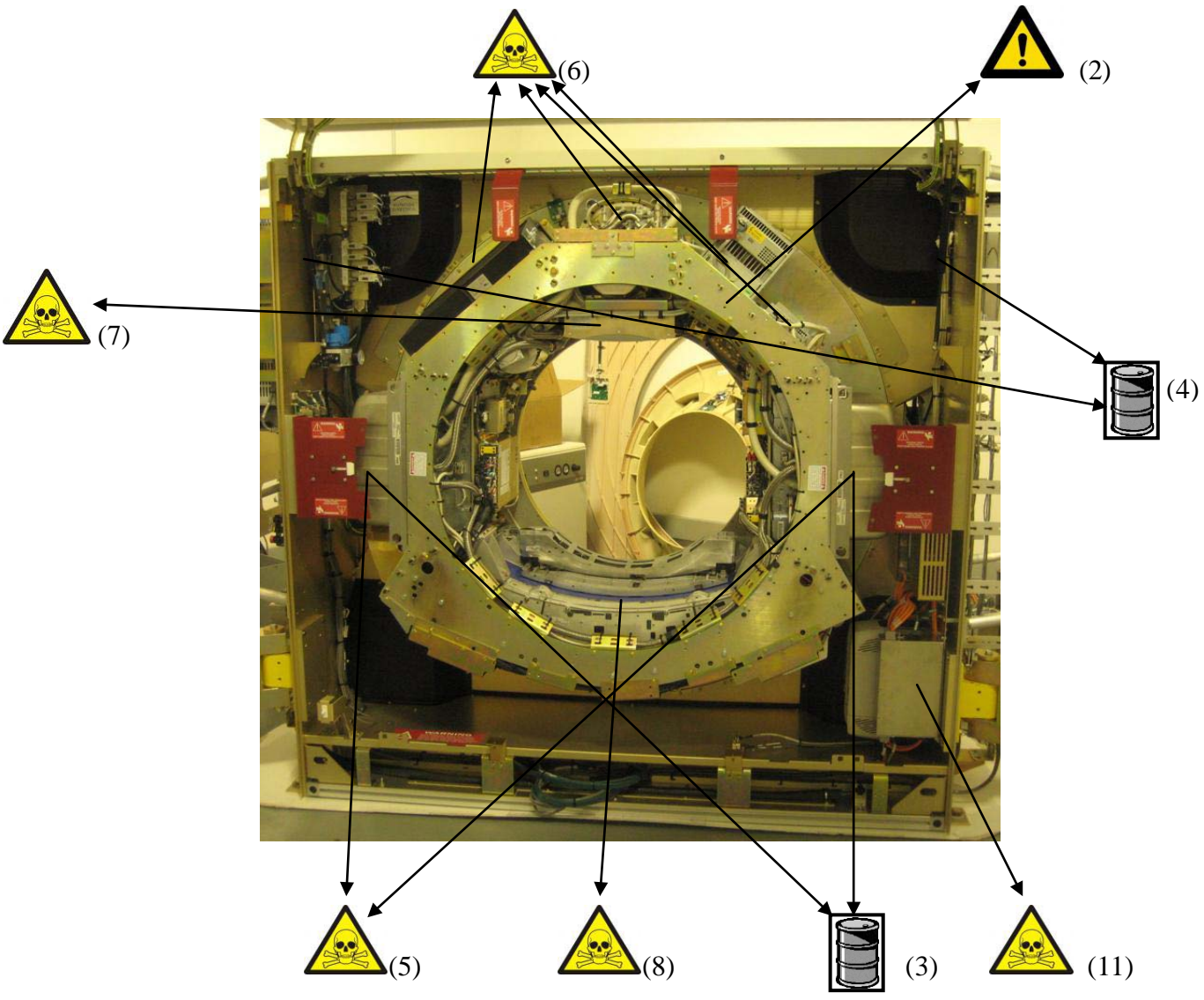


Figure 2: Front View of Gantry with Front Cover Raised

Title: Recycle Passport Brilliance iCT Gantry
Recycling Passport Number: 453567534481 Rev: B

Printed copies are uncontrolled

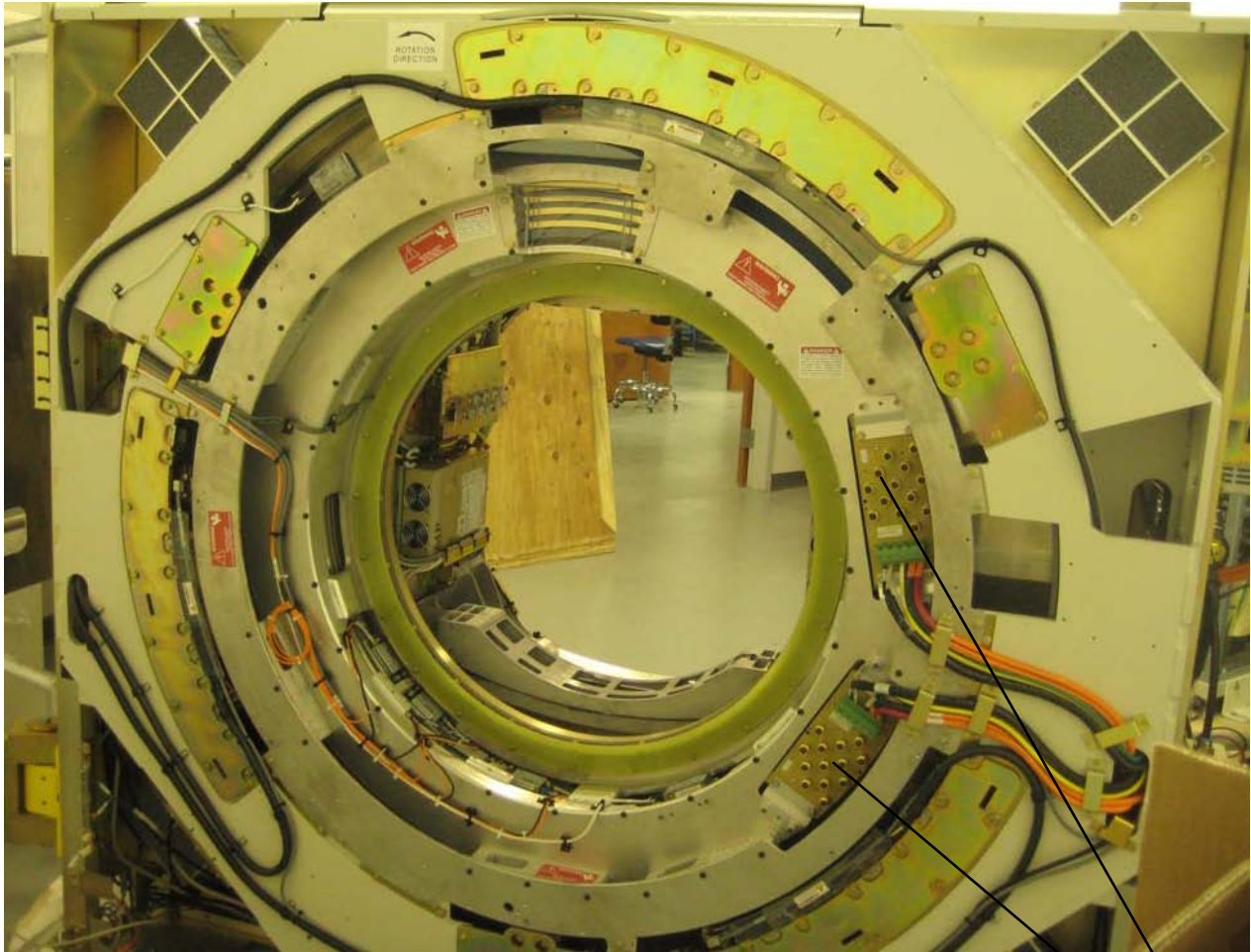


(9)

Figure 3: Back View of Gantry with Back Cover Removed

Title: Recycle Passport Brilliance iCT Gantry
Recycling Passport Number: 453567534481 Rev: B

Printed copies are uncontrolled

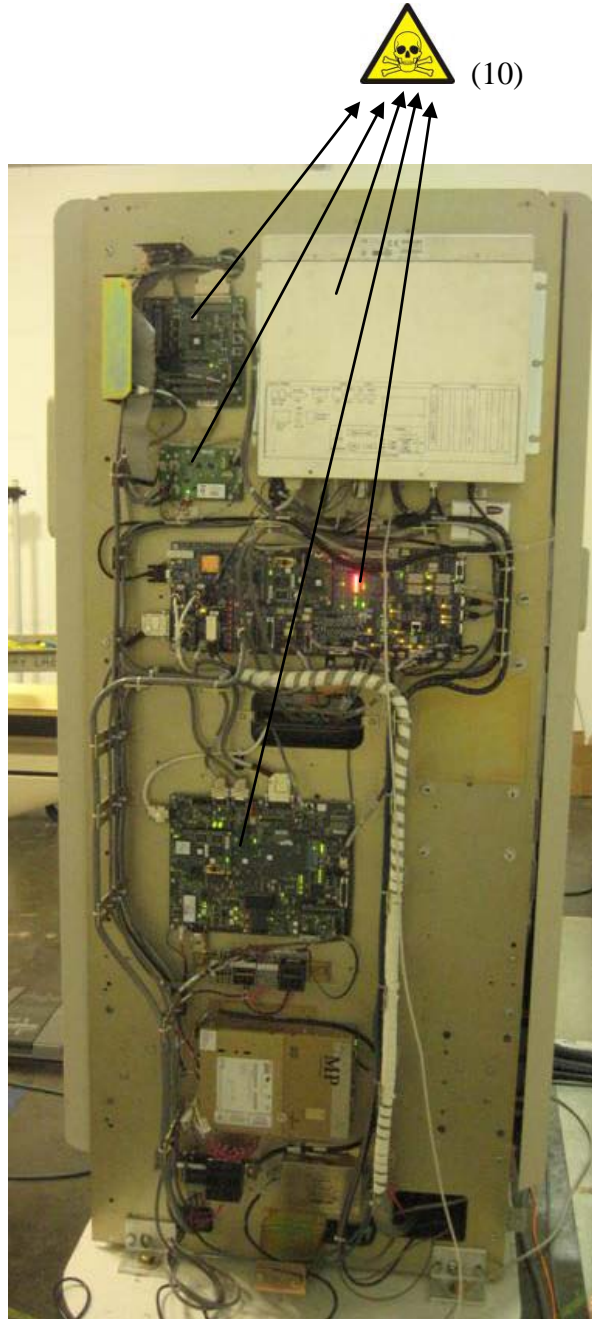

Title: Recycle Passport Brilliance iCT Gantry
Recycling Passport Number: 453567534481 Rev: B

Figure 4: Right Side View of Gantry with Column Cover Removed

Printed copies are uncontrolled

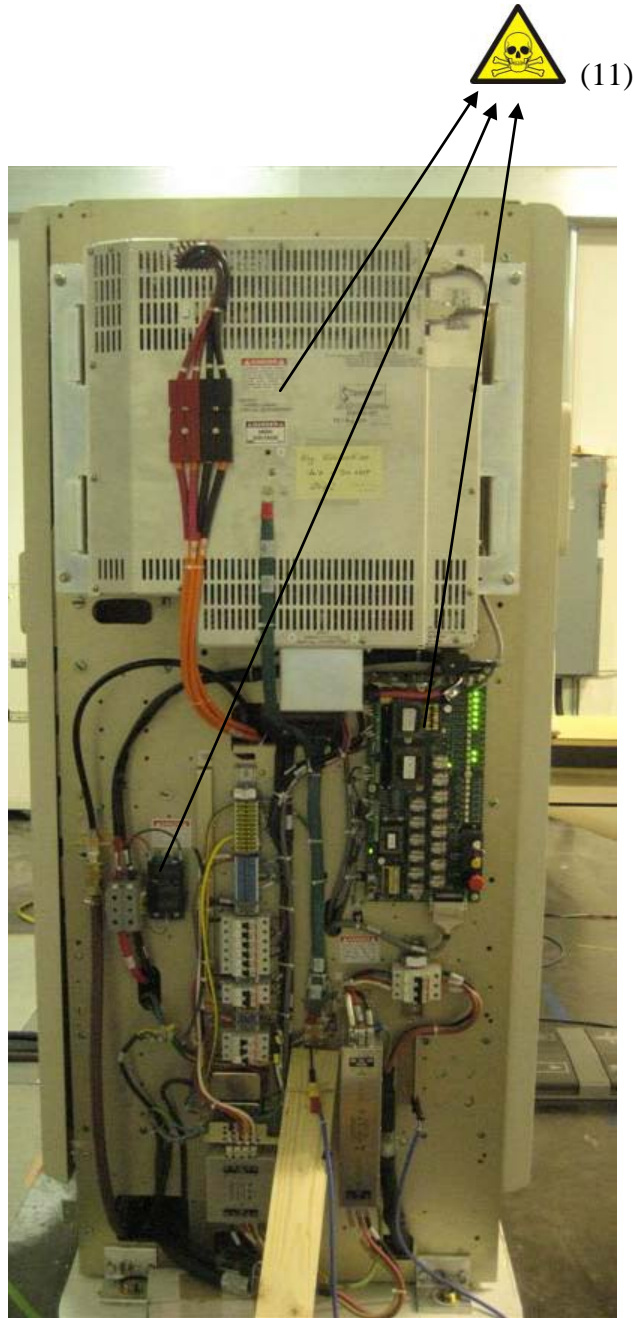


Figure 5: Left Side View of Gantry with Column Cover Removed

Title: Recycle Passport Brilliance iCT Gantry
Recycling Passport Number: 453567534481 Rev: B

PHILIPS HEALTHCARE	Product Recycling Passport	Page 7 of 7
-------------------------------	-----------------------------------	-------------

Recycling Passport Revision History

Rev	Description	Date
A	CO 2019848-Release Status A	24 February 2009
B	CO 2020759 - Added titles to images, removed superfluous footer and corrected Rev A description to be "CO 2019848"	9 June 2009

Printed copies are uncontrolled