
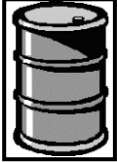





<b>Product name:</b>		<b>Ingenuity CT Gantry</b>
<b>Identification code(s)</b>		<b>459800029801</b>
<b>Total weight (in Kg)</b>		<b>1941</b>
<b>Producer/ Manufacturer</b>	Name company:	Philips Medical Systems
	Address:	595 Miner Road Highland Heights, Ohio
	Zip code:	44143
	Country:	United States of America


<b>Recycle Info</b>	<b>Items:</b>	<b>Location</b>
<b>Special attention</b> 	Spring loaded cover	Figure 1 (1)
	Actuators: Unstable when removed, possible rotation of frame	Figure 4 (2), Figure 5 (2)
	System Unbalanced when rotor components are removed causing rotation	Figure 3 (5)
<b>Fluids / Gases</b>		
	<b>Items:</b>	<b>Location</b>
	X-Ray Tube Hydrocarbon Oil	Figure 3 (2)
	Generators Hydrocarbon Oil	Figure 3 (6)
	Overflow Tank Hydrocarbon Oil	Figure 3 (7)
<b>Batteries</b>		
	<b>Type:</b>	<b>Location</b>
<b>To be Removed</b>		
<b>Hazardous</b>		
	<b>Substances:</b>	<b>Location</b>
	Printed Circuit Boards	Figure 1 (2), Figure 2 (1), Figure 4 (1), Figure 5 (1)
	Lead	Figure 3 (1),(3),(4)
<b>To be Removed</b>		

Title: Recycle Passport Ingenuity CT Gantry  
 Recycling Passport Number: 459800122671 Rev: A

Show locations of materials mentioned on the previous sheet. (Picture information)

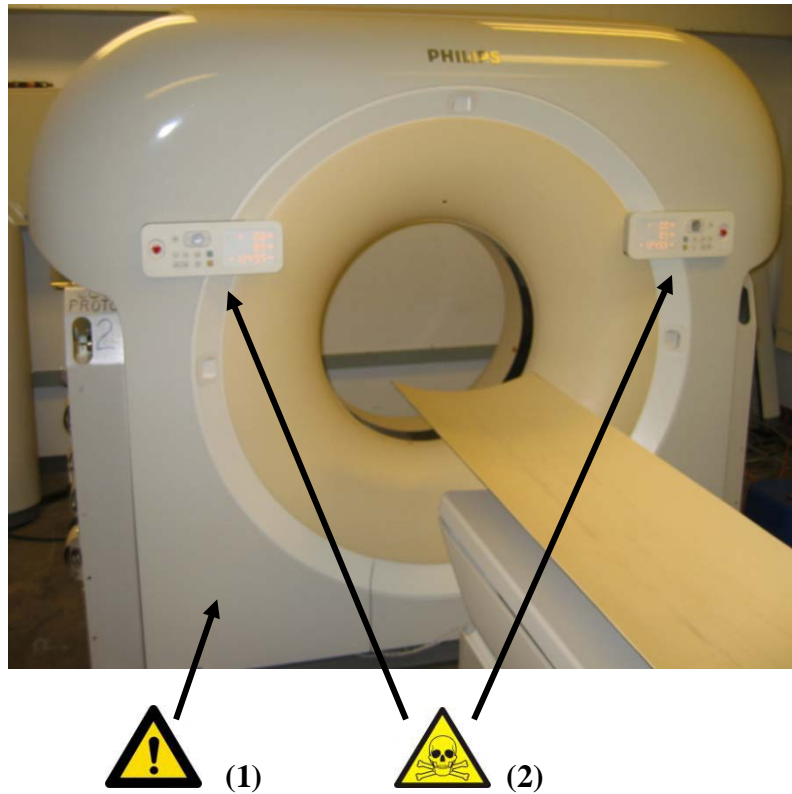


**Figure 1: Front View of Gantry with Covers on**

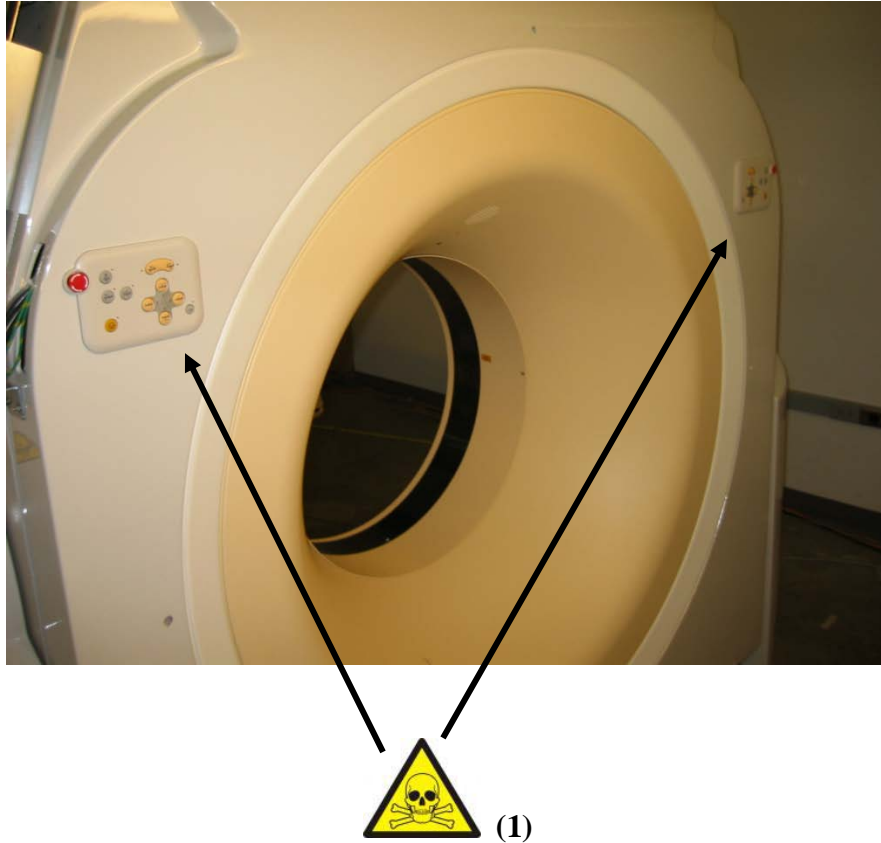


**Figure 2: Rear View of Gantry with Covers on**

Title: Recycle Passport Ingenuity CT Gantry  
Recycling Passport Number: 459800122671 Rev: A

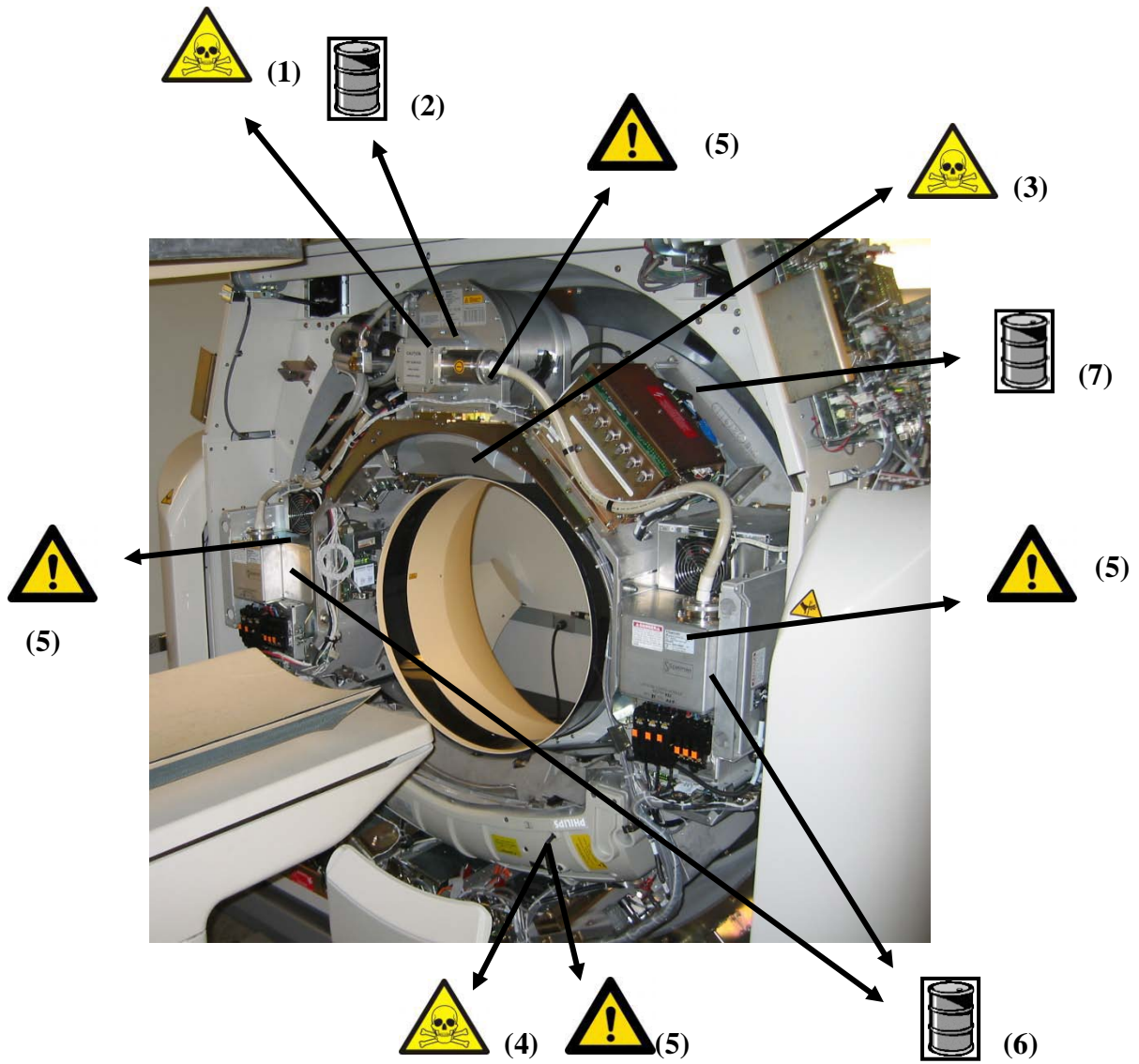


**Figure 3: Front View of Gantry with Front Cover Lifted**

Title: Recycle Passport Ingenuity CT Gantry  
Recycling Passport Number: 459800122671 Rev: A

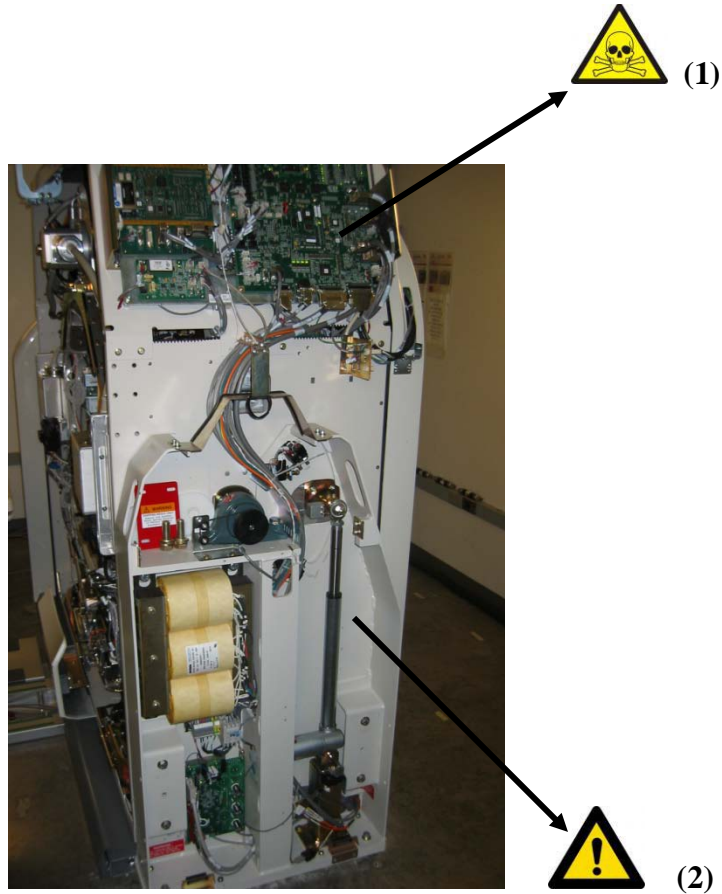

**Figure 4: Right Side View of Gantry with Cover Removed**

Title: Recycle Passport Ingenuity CT Gantry  
Recycling Passport Number: 459800122671 Rev: A

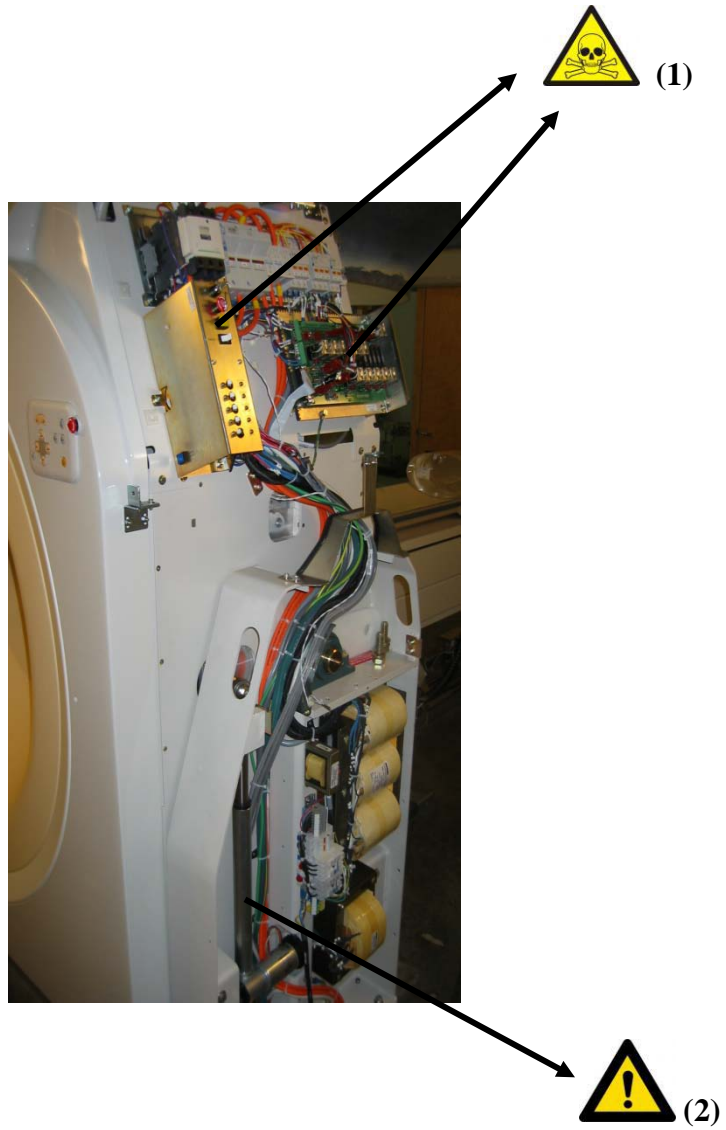

**Figure 5: Left Side View of Gantry with Cover Removed**

Title: Recycle Passport Ingenuity CT Gantry  
Recycling Passport Number: 459800122671 Rev: A