





Product Recycling Passport

Product name:	Precess Patient Monitor	
Identification code(s):	3160 (includes 3160DCU and 3160WPU)	
Total weight (in kg):	63 kg	
Producer /Manufacturer	Name company:	Invivo Corporation
	Address:	12601 Research Parkway Orlando, FL
	Zip code:	32826
	Country:	USA
	Electronic info:	Notice, Recycling Passport, 3160.doc

Recycle info	Item:	Location
	Special Attention	None
Fluids / Gases 	Type:	Location
	Not applicable	
Batteries  To be removed	Substances:	Location
	Batteries, Lithium-ion, Qty 4 (Invivo Part Nos. 9093 and 9064)	Figure 1 – Item 2 Figure 1 – Item 3 Figure 1 – Item 6 Figure 1 – Item 7
	Batteries, Lithium-polymer, Qty 2 (Invivo Part No. 9065)	Figure 6 – Item 1 Figure 7 – Item 1
	Batteries, Lithium , Qty 2	Figure 3 – Item 1 Figure 4 – Item 1

Product Recycling Passport

Hazardous	Item:	Location
 To be removed	Lead (Pb) – Printed circuit board solder	Figure 1 – Item 1 Figure 5 – Item 1 Figure 6 – Item 3 Figure 7 – Item 3 Figure 8 – Item 1
	Lead (Pb) – Patient cables	Figure 6 – Item 2 Figure 7 – Item 2 Figure 10 – Item 1
	Lead (Pb) – Fiber optic cables	Figure 2 – Item 2 Figure 9 – Item 1
	Lead (Pb) – Display	Figure 1 – Item 4
	Hexavalent Chromium (Cr ⁶⁺) – Metal enclosures	Figure 1 – Item 5 Figure 2 – Item 1 Figure 5 – Item 2
	Cadmium (Cd) – Patient cables	Figure 6 – Item 2 Figure 7 – Item 2 Figure 10 – Item 1

Product Recycling Passport

Locations as mentioned in the Passport (Pictures information)



Figure 1: 3160 Precess Patient Monitor

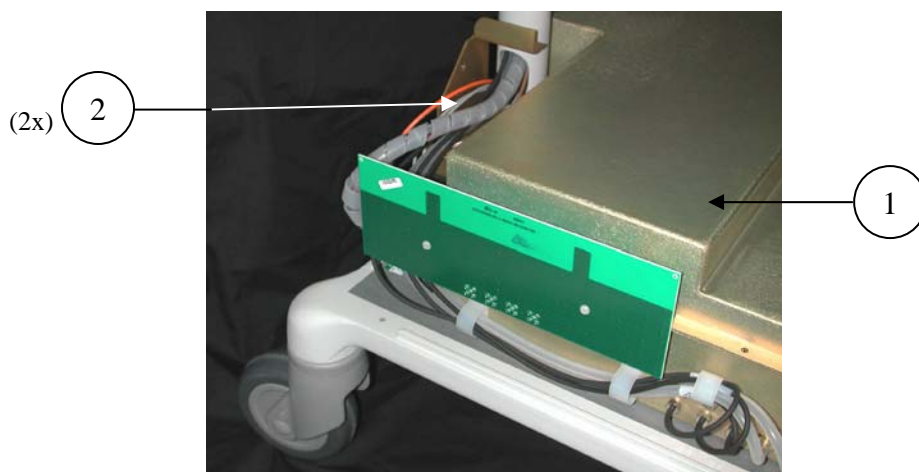


Figure 2: 3160 Precess Patient Monitor (WPU, Internal View)

Product Recycling Passport

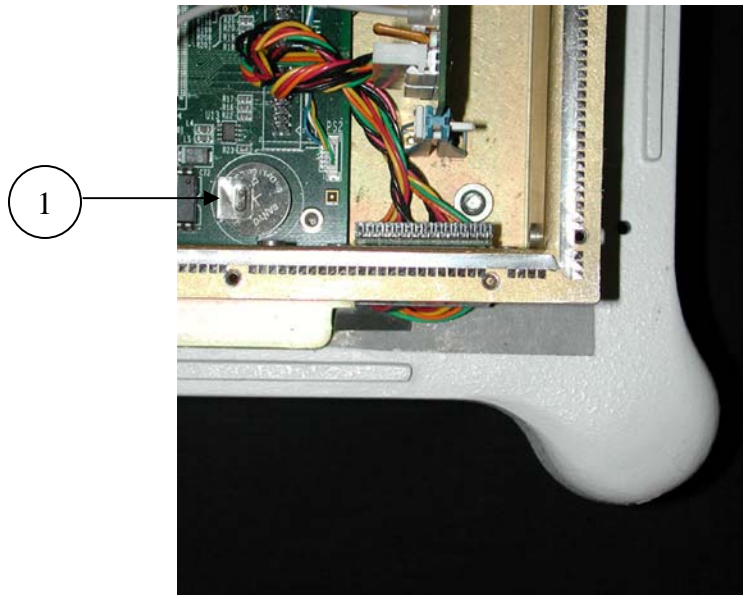


Figure 3: WPU, View of Power Supply Printed Circuit Board



Figure 4: DCU, View of Power Supply Printed Circuit Board

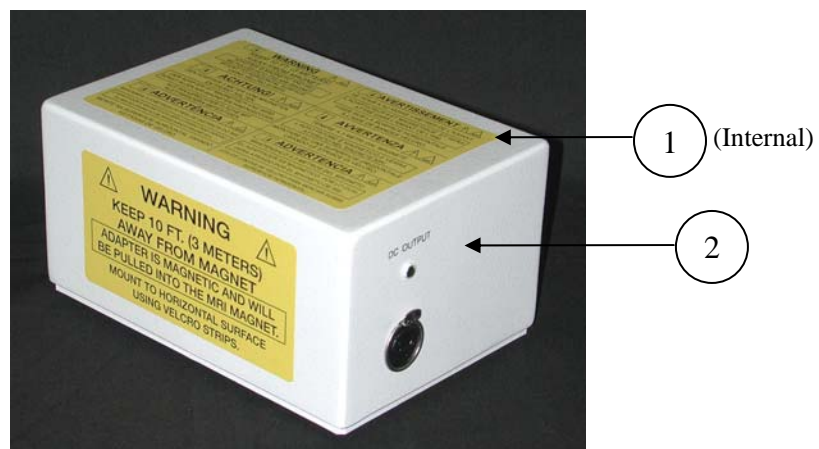
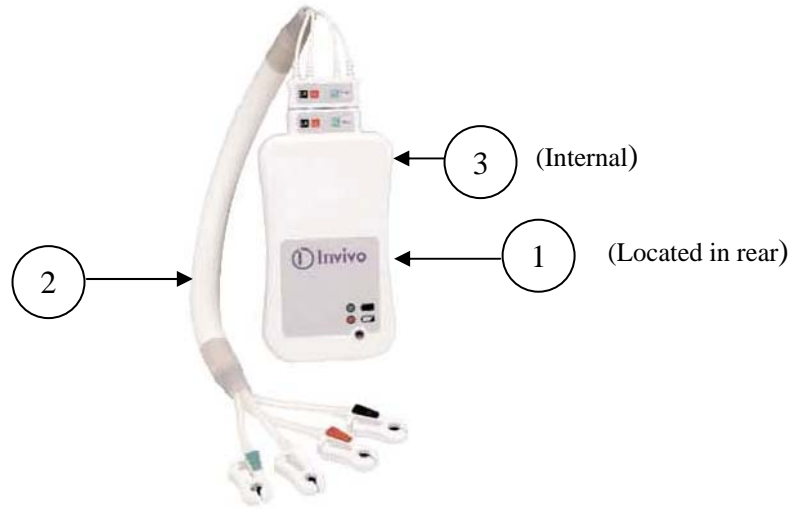
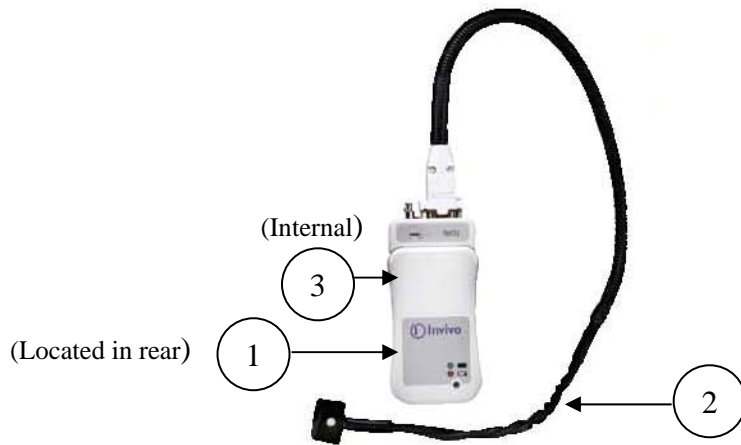


Figure 5: External Power Adapter (Invivo Part No. AS201)

Product Recycling Passport



**Figure 6: WECG Module with ECG Leads
(Invivo Part Nos. 9312 and 9224)**



**Figure 7: WSPO2 Module with SPO2 Leads
(Invivo Part Nos. 9311 and 9399B)**

Product Recycling Passport



Figure 8: Battery Charger (Invivo Part No. 9023)



Figure 9: Fiber-Optic Surface Temperature Probe (Invivo Part No. 9320)



Figure 10: NIBP Kit (Invivo Part No. 94014)