

## Product (*IntelliBridge Hubs*) Recycling Passport

Notice                    **UNCONTROLLED IN PRINTED FORM UNLESS STAMPED IN RED**

---

**Document  
Ownership**            The document owner is listed below:

<b>Function</b>	<b>Title</b>
<i>PMD-A Supply Chain Engineering Department</i>	<i>Manufacturing Engineer</i>




---

**Revision  
Summary**            The current document revision is described below:

<b>Rev.</b>	<b>Change Order #</b>	<b>Change Summary</b>
<i>A</i>	<i>See eDMR</i>	<i>As Issued</i>

## Product (*IntelliBridge Hubs*) Recycling Passport

<b>Product name:</b>	<i>Enter name of product as it appears on the product.</i> 453564063601, 453564063611 IntelliBridge EC40/80 Hub	
<b>Reference Identifier:</b>	<i>This should be the model number as it appears on the product.</i> 865056	
<b>Total weight (in Kg):</b>	1.4Kg	
<b>Producer /Manufacturer</b>	Name company:	Philips Healthcare
	Address:	3000 Minuteman Road
	Zip code:	Andover, MA 01810
	Country:	USA

Recycle info	Substances:	Location
<b>Special attention</b> 	<i>Indicate any part of the product that could pose a hazard to the recycler and treatment facility during disassembly such as capacitors with residual energy or spring loaded valves and show location</i>	N/A
<b>Batteries</b>  <b>To be removed</b>	Type:	Location
	<i>Indicate all batteries in the product via battery chemistry and show location</i>	N/A
<b>Hazardous</b>  <b>To be removed</b>	Item	Location
	<i>Indicate location of printed circuit boards and mercury in back lights as well as any cadmium, flame retardants (PBB and PBB) and hexavalent chromium in the product and show location</i>	
	General disassembly	See figure 1
	453564063621 System PCA	See figure 2
	453564063661 Option PCA	See figure 3
	453564063681 LED PCA	See figure 4
453564081181 Power Supply	See figure 5	

### Additional Information as needed

Material content (global)	Weight in Kg	Materials content (continue)	Weight in Kg
453564063621 System PCA	0.160		
453564063661 Option PCA	0.095		
453564063681 LED PCA	0.025		
453564081181 Power Supply	0.170		

## Product (*IntelliBridge Hubs*) Recycling Passport

Locations as mentioned in the Passport (Pictures information)

**Figure #1: Remove Plastic Cosmetic Cover**



Push cover upward



Plastic Cosmetic cover



Remove screws (10)

## Product (*IntelliBridge Hubs*) Recycling Passport

Figure #2: Removing the Option PCA



Remove connector

Remove screws (4)



453564063661 Option PCA

## Product (*IntelliBridge Hubs*) Recycling Passport


Figure #3: Removing the System PCA



Remove screws (6)



Remove cable  
Push outward

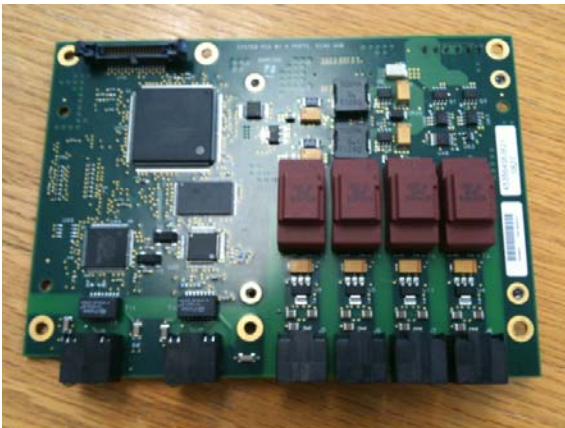

453564063621 System PCA

## Product (*IntelliBridge Hubs*) Recycling Passport

Figure #4: Removing the LED PCA



Remove screws (4)



453564063681 LED PCA

## Product (*IntelliBridge Hubs*) Recycling Passport

Figure #5: Removing the Power Supply



Remove screws (4)



453564081181 Power Supply