
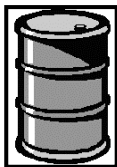



PHILIPS	Product Recycling Passport LC0593 – WEEE Passport, FR3 Rechargeable Battery	Page 1 of 4
----------------	--	-------------

Product name:	Assy, Battery, Rechargeable Use, FR3	
Identification code(s)	453564071819; 989803150241	
Total weight (in Kg)	0.63 kg -0.65 kg	
Producer/ Manufacturer	Name company:	Philips Medical Systems
	Address:	22100 Bothell Everett Highway
	Zip code:	Bothell, WA 98021
	Country:	United States

Recycle Info	Items:	Location
Special attention 	Remove welded cover from Battery Enclosure	Figure 1 & 2
	Contains circuit board.	Figure 4, 5 & 6

	Items:	Location
 Fluids / Gases	N/A	

	Type:	Location
Batteries  To be Removed	Li Ion Cell, 18650 size AA, quantity: 6	Figure 3


	Substances:	Location
Hazardous  To be Removed	N/A	


FIGURE 1**RECHARGEABLE BATTERY, TOP VIEW****FIGURE 2****USING SIDE CUTTER TOOL, SNIP EACH CORNER AT WELD LINE (SEAM) – FOUR PLACES**

FIGURE 3



PRY COVER OPEN TO EXPOSE CELLS.
REMOVE FROM PLASTIC CASE

FIGURE 4

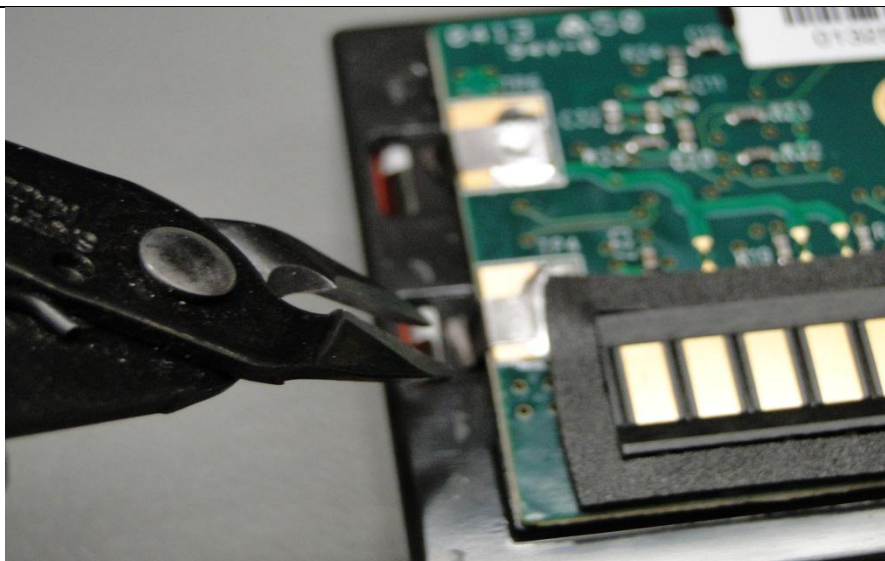

CUT NICKEL CONTACT STRIPS TO DISCONNECT CIRCUIT BOARD FROM
PLASTIC BRACKET, FOUR PLACES
(SEE FIGURE 5)

FIGURE 5



NICKEL CONTACT LOCATIONS

FIGURE 6



SCREW

**USING #10 TORX DRIVER, REMOVE SCREW TO FREE CIRCUIT BOARD
 FROM PLASTIC BRACKET**

