
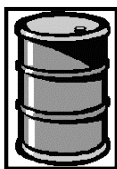




Product name:	Q-CPR Bundle Box Assembly	
Identification code(s)	989803149941	
Total weight (in Kg)	1.2 kg	
Producer/ Manufacturer	Name company:	Philips Medical Systems
	Address:	22100 Bothell Everett Highway
	Zip code:	Bothell, WA 98021
	Country:	United States

Recycle Info	Items:	Location
Special attention 	Q-CPR Meter Assembly Contains printed circuit assembly (PCA) that is not readily accessible to end-user.	Figures 1 & 2 Figures 2, 3 & 4
 Fluids / Gases	N/A	
Batteries  To be Removed	N/A	
Hazardous  To be Removed	LCD (display) contains mercury (Hg)	Figure 1

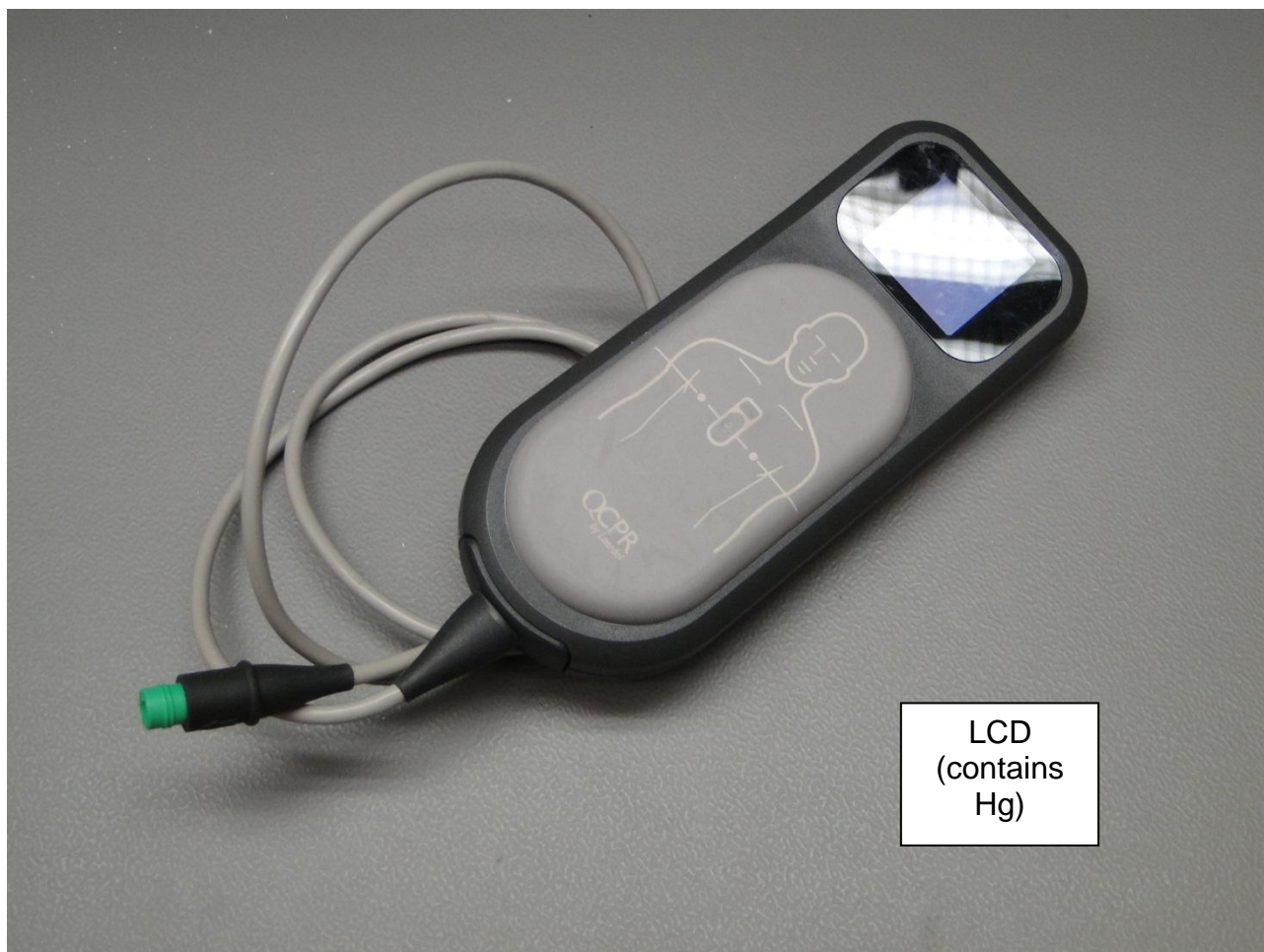
**FIGURE 1**

FIGURE 2



Q-CPR Meter, BACK VIEW

2 SCREWS

**USING STRAIGHT-BLADE SCREWDRIVER, REMOVE (2) SCREW VIEW INTERIOR
CAVITY OF Q-CPR METER**



FIGURE 3

**NO BATTERY CELL IS INSTALLED IN THIS DEVICE
AS CONFIGURED TO FUNCTION WITH PHILIPS FR3 AED.**

THERE IS A PRINTED CIRCUIT BOARD ASSEMBLY (PCA) LOCATED UNDER THE BATTERY CAVITY.

**NOTE: THIS DEVICE IS AN ULTRASONICALLY WELDED ENCLOSURE.
THERE IS NO READILY- AVAILABLE ACCESS TO THE PCA INSIDE THIS DEVICE.**



IF REQUIRED BY END-USER, ACCESS TO THE PCA AND THE LCD
INSIDE THIS UNIT CAN BE GAINED BY CUTTING THE DEVICE OPEN ALONG
THE WELDED SEAM LINE AS INDICATED IN FIGURE 3.

FIGURE 4

IF REQUIRED, WELDED SEAM WILL NEED TO BE MACHINE-CUT
TO ALLOW ACCESS TO INTERNAL PCA.