

# SureSigns VS3

Vital signs monitor



## Innovative vital signs

The SureSigns VS3 vital signs monitor provides caregivers with the information they require in an affordable and easy-to-use package. Built on our strong heritage and industry leadership in patient monitoring, all Philips SureSigns monitors blend vital measurement and monitoring capabilities with simplicity, convenience, and reliability.





The programmable blood pressure feature allows you to create up to five different programs to monitor NBP at particular intervals based on your own unit-specific protocols.

### **Better patient care**

Because the SureSigns VS3 user interface is intuitive and easy to use, you can focus on what matters most—patient safety and comfort. The full-color screen contains large, easy-to-read numerics, a scrollable list of stored patient records, and an optional  ${\rm SpO}_2$  waveform display. The innovative pop-up screens, fixed keys, and icon-based menus provide quick access to frequently used features.

Available SureSigns VS3 measurements include blood pressure, pulse oximetry, pulse rate, and temperature. Designed for a wide variety of care environments, the SureSigns VS3 is suitable for the general care ward, emergency department triage, out-patient surgery, urgent care centers, or anywhere else you need basic vital signs monitoring.

### Superior quality and support


The SureSigns VS3 monitor comes with the quality you expect from Philips and a design that is easy to support. Every monitor includes a lithium ion battery and AC power cord clip to enhance battery operation and charging performance. Software upgrades, via the USB port, are remarkably quick and easy.

Replacement parts that are compatible with the SureSigns VM monitors help to streamline parts inventory management, and the sleek modular design simplifies monitor repair. As always, Philips monitors come with a world-class warranty and 24-hour\* telephone support from our global response centers.



### monitor.

### Easy to use. Easy to support.

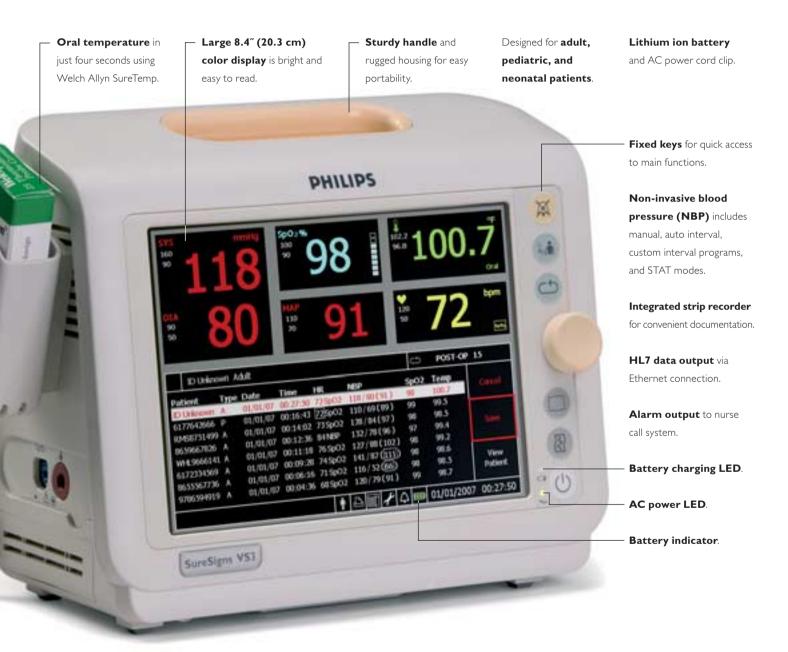

Patient ID entry via **optional** 

bar code scanner.

**USB port** for easy software upgrades.

Philips SpO, technology

is compatible with a wide variety of adult, pediatric, and neonatal sensors from Philips, Nellcor,\*\* and Masimo.\*\* SpO, waveform display.

Stores up to 400 patient records.

Product specifications			
Width, height, depth	W 26 cm (10.2"), H 22 cm (8.6"), D 14.5 cm (5.7")		
Weight	3.1 kg (6.9 lbs), excluding optional recorder		
Screen	20.3 cm (8.4") SVGA TFT-AM LCD		
NBP	Oscillometric using stepwise deflation pressure Adult/pediatric measurement range - Systolic: 30–255 mmHg (4–34 kPa) - Diastolic: 15–220 mmHg (2–29 kPa) - Mean: 20–235 mmHg (2.7–31.3 kPa) Neonatal measurement range - Systolic: 30–135 mmHg (4–18 kPa) - Diastolic: 15–110 mmHg (2–14.7 kPa) - Mean: 20–125 mmHg (2.7–16.7 kPa)  • Blood pressure accuracy: ±5 mmHg (.70 kPa) • Pulse rate accuracy: ±2 bpm or ±2% (whichever is greater)		
SpO <sub>2</sub>	<ul><li>Meets EN ISO 9919 standard</li><li>Measurement range: 0–100%</li><li>Accuracy: depends on sensor</li></ul>		
Temperature	<ul> <li>Modes: Predictive and Monitored</li> <li>Probe sites: oral, rectal, or axillary</li> <li>Range: 26.7–43.3°C (80–110°F)</li> <li>Accuracy: ±0.2°C (±0.4°F)</li> <li>Meets EN12470-3 standard</li> </ul>		
Battery	<ul> <li>Lithium ion, 11.1 V</li> <li>Operating time: 4 hours with NBP every 15 minutes.</li> </ul>		
Data output	HL7 format, via Ethernet port		
Patient type	Adult, pediatric, neonatal		

Standard accessories	To re-order
Adult antimicrobial NBP cuff	M4555A
NBP interconnect tubing	M1599B
Adult finger SpO <sub>2</sub> sensor	M1196A
Oral/axillary temperature probe kit	989803143381
Temperature probe covers (one box)	M4823A
Recorder paper (5 rolls)	989803136891
A Little III	Blair

Additional accessories are available.	Please contact your Philips representative for
more information.	

Also available	To order
Roll stand with basket	989803144001
Wall mount	989803144011
Bar code scanner kit	989803147821



**Supplies for the SureSigns VS3** include a wide selection of disposable and reusable SpO<sub>2</sub> sensors and various types and sizes of blood pressure cuffs—all designed for maximum patient comfort, performance, and reliability. A starter kit of Philips supplies is included with every order.


**Quick unpack and install**. Each monitor ships with a set-up guide, instructions for use, quick reference cards, a user training DVD, and a service manual on CD.



\* Only available in US.

\*\* Adapter cable required.

buy online\*
http://shop.medical.philips.com

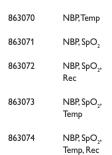

NBP only

863069

© 2006 Koninklijke Philips Electronics N.V.

All rights are reserved. Philips Medical Systems Nederland B.V. reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

Nellcor is a registered trademark of Nellcor Puritan Bennett, Inc. Masimo is a registered trademark of Masimo Corporation. SureTemp is a registered trademark of Welch Allyn.

Philips Medical Systems is part of Royal Philips Electronics

www.medical.philips.com medical@philips.com fax: +31 40 27 64 887

Printed in The Netherlands 4522 962 20051/862 \* DEC 2006 Philips Medical Systems Global Information Center P.O. Box 1286 5602 BG Eindhoven The Netherlands

