
Product name: iE33		iE33 M	edical Diagnostic Ultrasound Equipment	
Identification code(s)		8500-0082-01		
Total weight (in Kg)		163 Kg		
	N			
Producer /	Name company:		Philips Ultrasound	
Manufacturer	Address:		22100 Bothell Everett Highway	
	Zip code:		98041-3003	
	Country:		USA	


Recycle Info	Items:	Location
Special attention	Gas Spring under pressure. See Figure 1 [1]	
Fluids / Gases	Items:	Location
	Oil and Nitrogen Gas: (2X) in Gas Spring Tubes in Control Panel Articulation Assy Do not drain oil with device under pressure.	Figure 1 [1]
	Oil and Nitrogen Gas: in Gas Spring Tube in Arm Articulation Assy Do not drain oil with device under pressure.	Figure 1 [2]
Batteries	Туре:	Location
Batteries	Li Battery: (2X) on CPU PCB	Figure 1 [3]
	Battery Lead Acid (Pb): in Power Supply Assy	Figure 1 [4]
To be Removed		
Hazardous	Substances:	Location
\wedge	Lead (Pb) – contains lead in PCB's, solder, wiring, and cables	N/A
	Mercury (Hg) – Video monitor lamps	Figure 2 [5]
No.	Mercury (Hg) – (2X) Touch screen lamps	Figure 2 [6]
To be Removed	Note: Transducers (Probes) contain lead (Pb) on Flex circuits / PCB's in the connector and probe areas	N/A

Title: iE33 Recycling Passport Recycling Passport Number: 8500-0082-00047 Rev: B

Figure 1



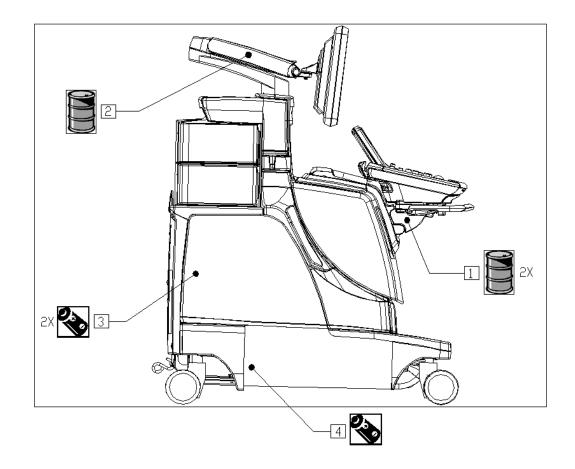

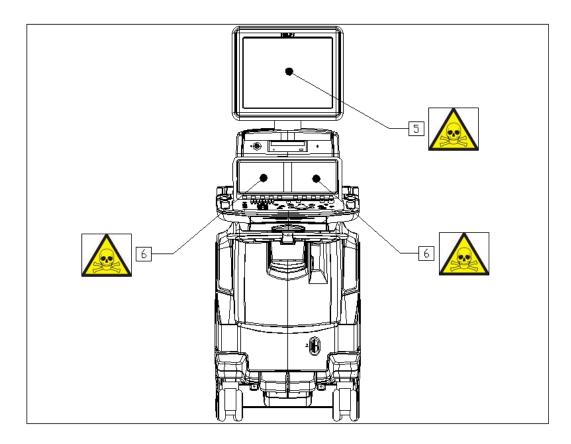

Philips Medical Systems/ Product Recycling Passport

Recycling Passport Revision History

Rev	Description	Date
А	First Passport	July 31, 2006
В	Revision for Product Weight Correction	September 14, 2009

Title: iE33 Recycling Passport Recycling Passport Number: 8500-0082-00047 Rev: B