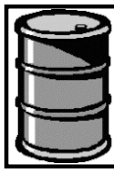


PHILIPS HEALTHCARE	Product Recycling Passport	Page 1 of 5
-------------------------------	-----------------------------------	-------------

Product name:	Brilliance iCT Extended Patient Support	
Identification code(s)	459800047821	
Total weight (in Kg)	408	
Producer/ Manufacturer	Name company:	Philips Medical Systems
	Address:	595 Miner Road Highland Heights, Ohio
	Zip code:	44143
	Country:	United States of America

Recycle Info	Items:	Location
Special attention 	Removal of Motor: Couch will drop when motor is removed.	Figure 1 (1)
	Green Safety Block: When removed, couch may collapse.	Figure 1 (3)
	Giant Spring: Potential energy hazard.	Figure 1 (4)
Fluids / Gases 	Items:	Location
Batteries 	Type:	Location
To be Removed		
Hazardous 	Substances:	Location
	Printed Circuit Boards	Figure 1 & 2 (2)
To be Removed		

Title: Recycle Passport Brilliance Extended Patient Support
 Recycling Passport Number: 459800081021 Rev: A

Show locations of materials mentioned on the previous sheet. (Picture information)

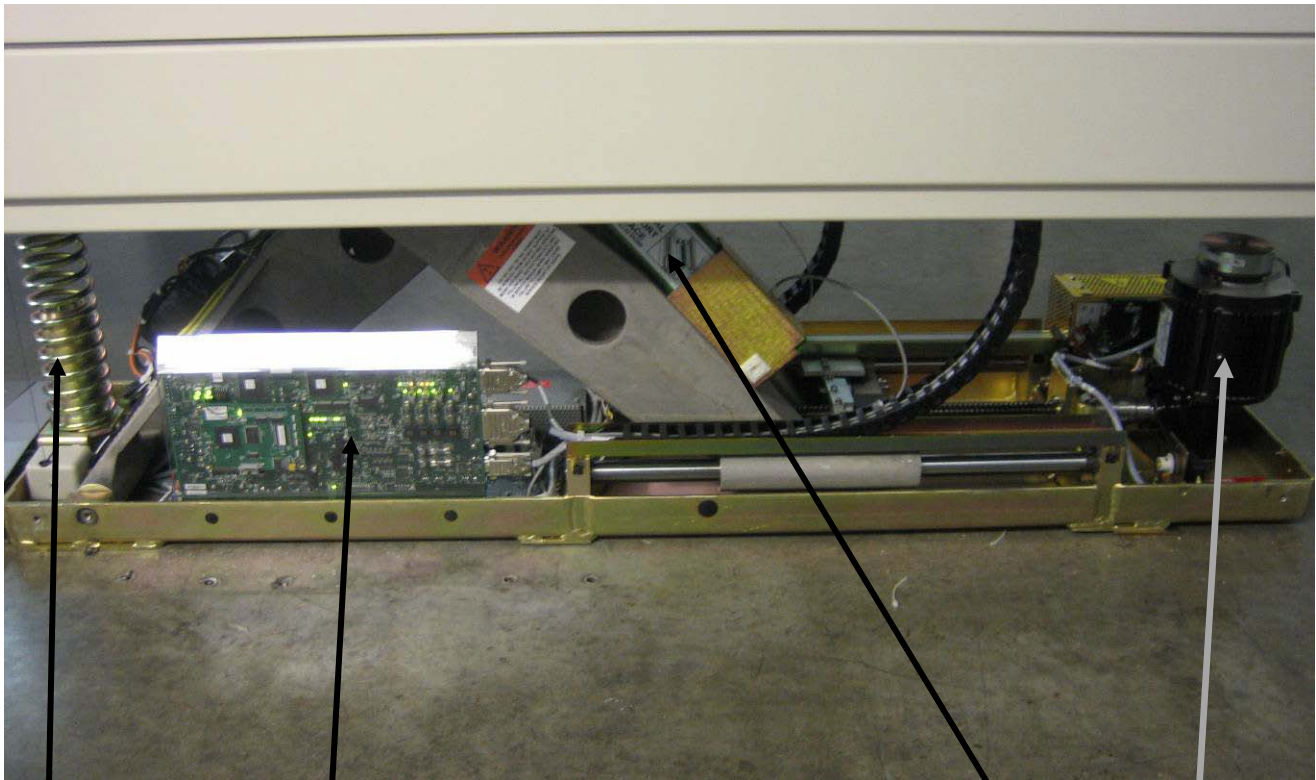


Figure 1: Right Side of Couch



(4)



(2)



(3)



(1)

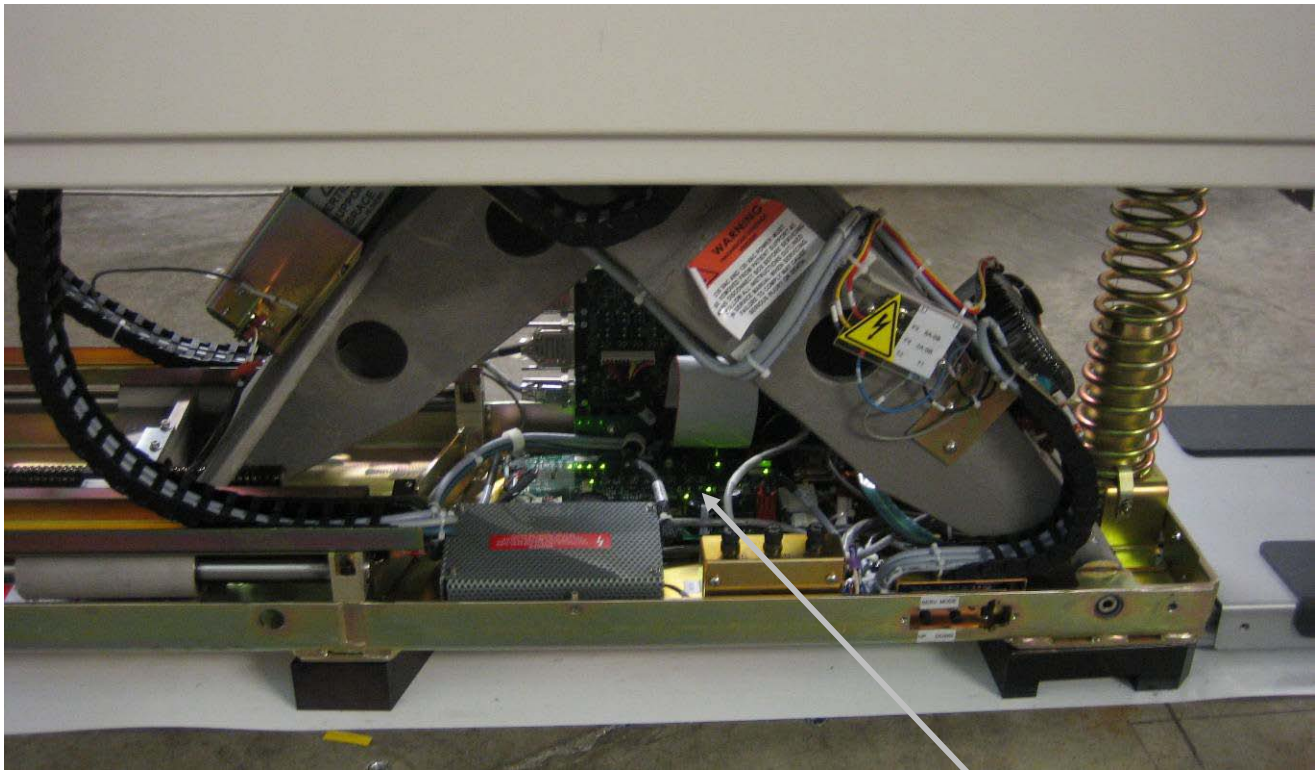


Figure 2: Left Side of Couch



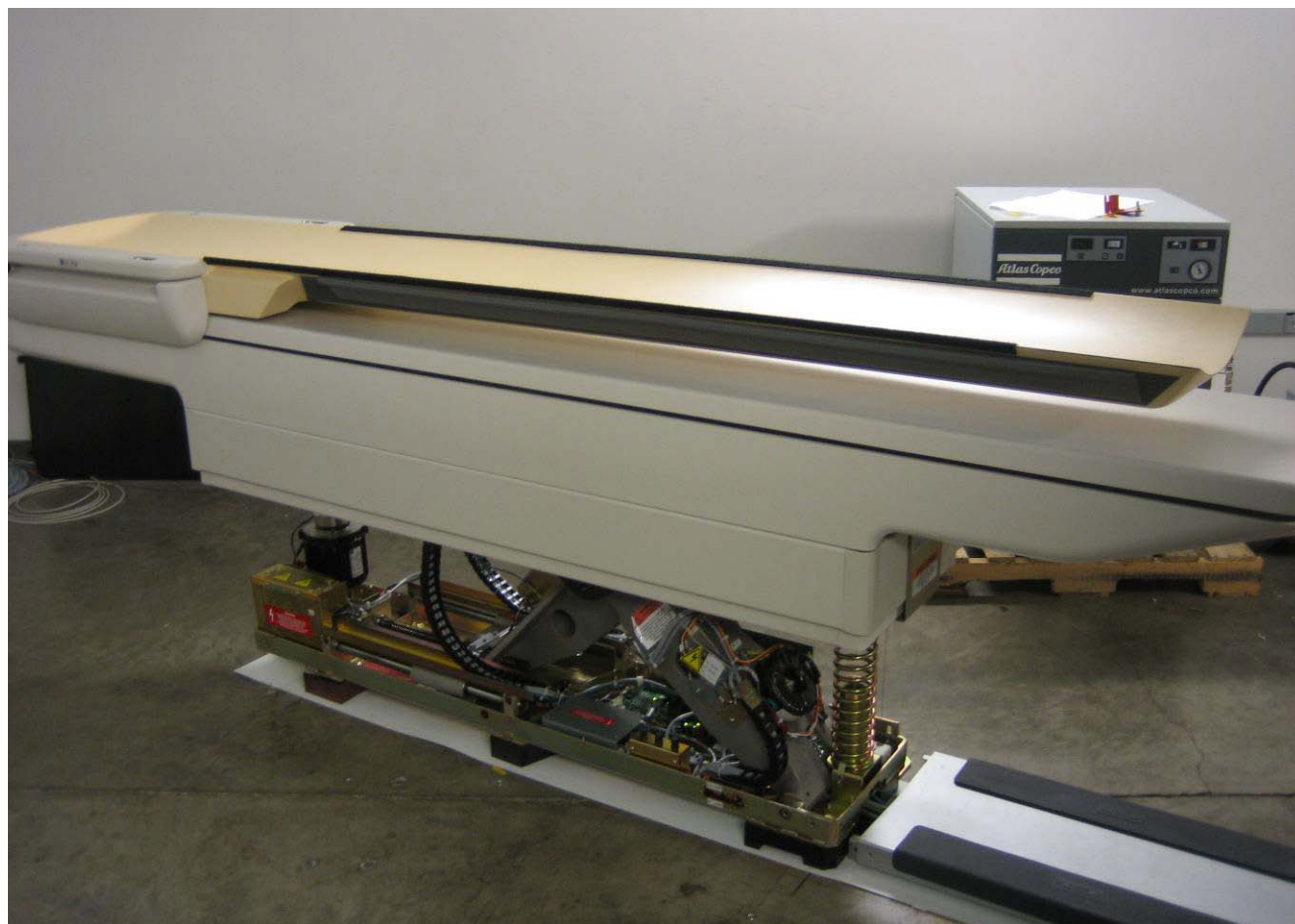


Figure 3: Full view of Patient Support

Title: Recycle Passport Brilliance Extended Patient Support
Recycling Passport Number: 459800081021 Rev: A

PHILIPS HEALTHCARE	Product Recycling Passport	Page 5 of 5
-------------------------------	-----------------------------------	-------------

Recycling Passport Revision History

Rev	Description	Date
A	E-007263-Release Status A	18 August 2010

Printed copies are uncontrolled